

User Guide

Smart connected toothbrush
Kolibree App iOS, Android

EN

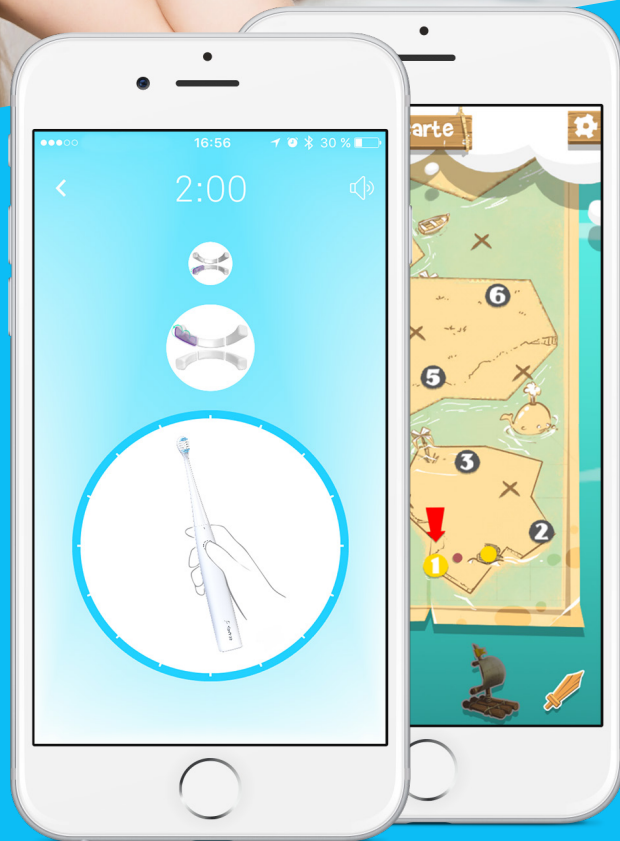



Table of contents

1. Getting started.....	P.4
2. Connecting the kolibree toothbrush.....	P.7
3. Brushing teeth.....	P.12
4. Kolibree account and users	P.15
5. Kolibree App.....	P.22
6. Brushing data.....	P.30
7. Go! Pirate Game.....	P.31
8. Kolibree Coach Activity	P.38
9. Kolibree toothbrush maintenance.....	P.41
10. Technical information.....	P.46
11. Important precautions.....	P.48
12. Legal information.....	P.50
13. Sending comments.....	P.53

1. Getting started


Description

- 1 - Kolibree interchangeable brush
- 2 - Button
- 3 - Kolibree toothbrush
- 4 - Charging base.

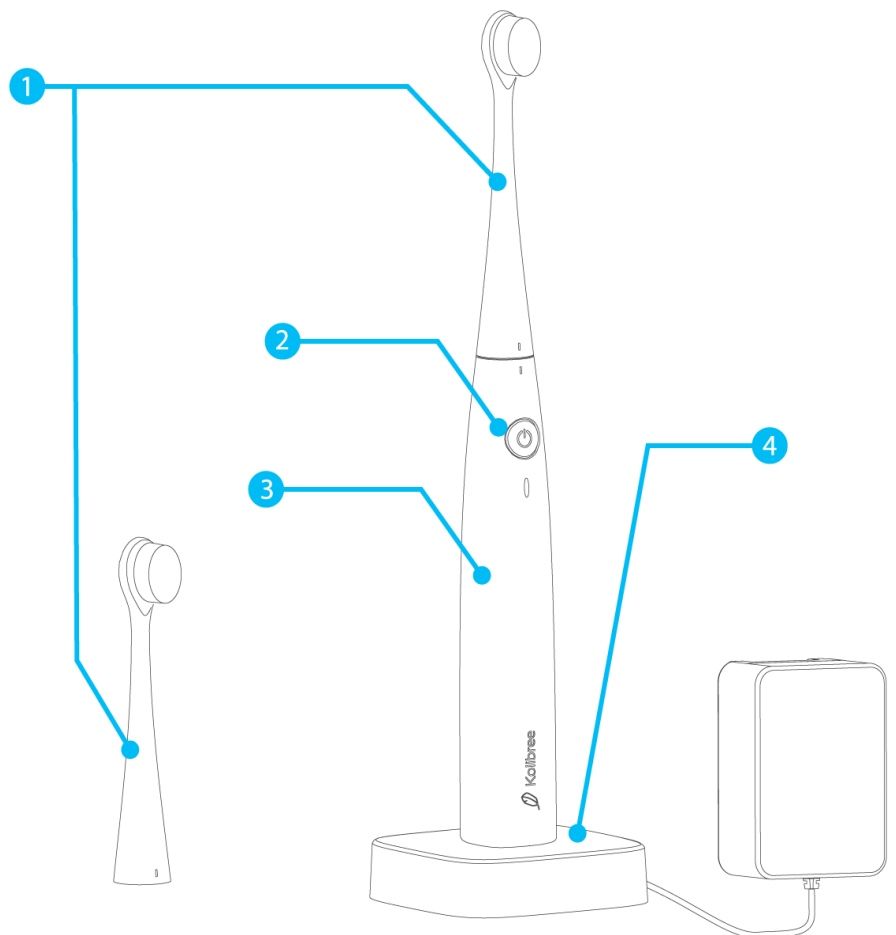

Charging the toothbrush


Your Kolibree toothbrush comes with an induction charger. To charge the battery, connect it to an outlet and place the toothbrush on the base. It charges upon contact.

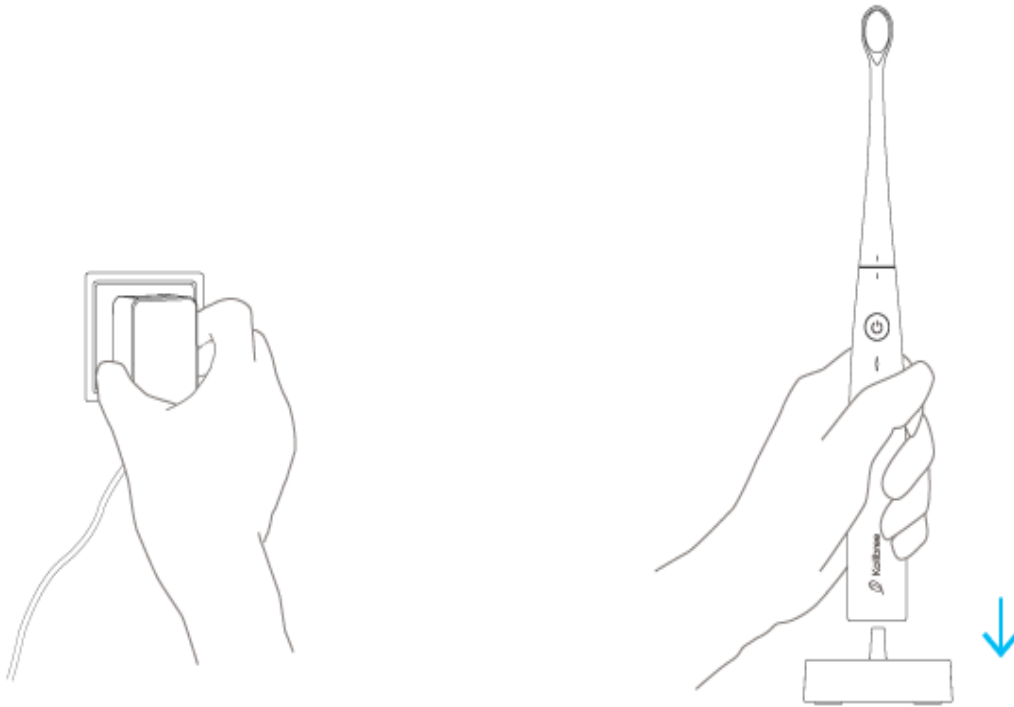

On/Off button

When turning on your Kolibree toothbrush for the first time, press the button once. The indicator light will briefly flash green for a few seconds.


Then press the button once to turn on the toothbrush, and the toothbrush will start to vibrate. To stop the toothbrush's vibration, press the button once again.

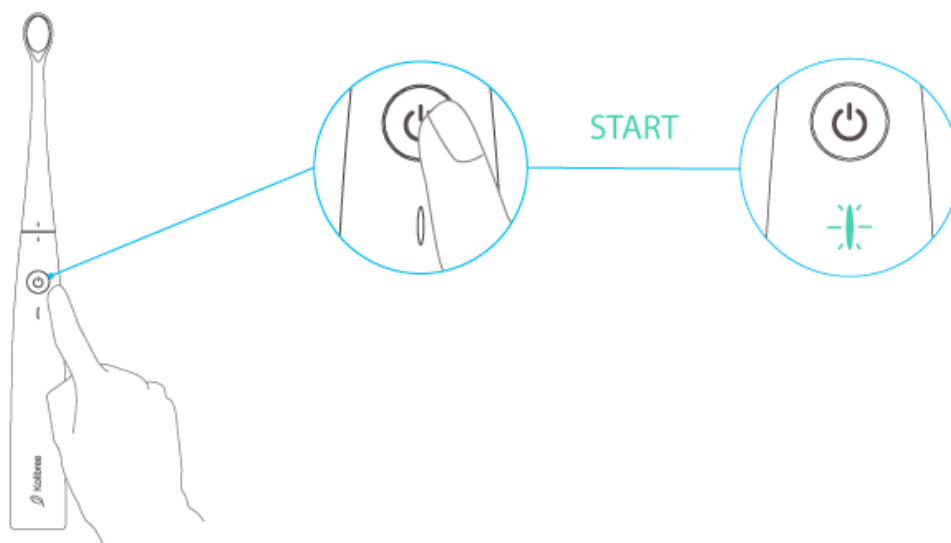

Turning off the toothbrush



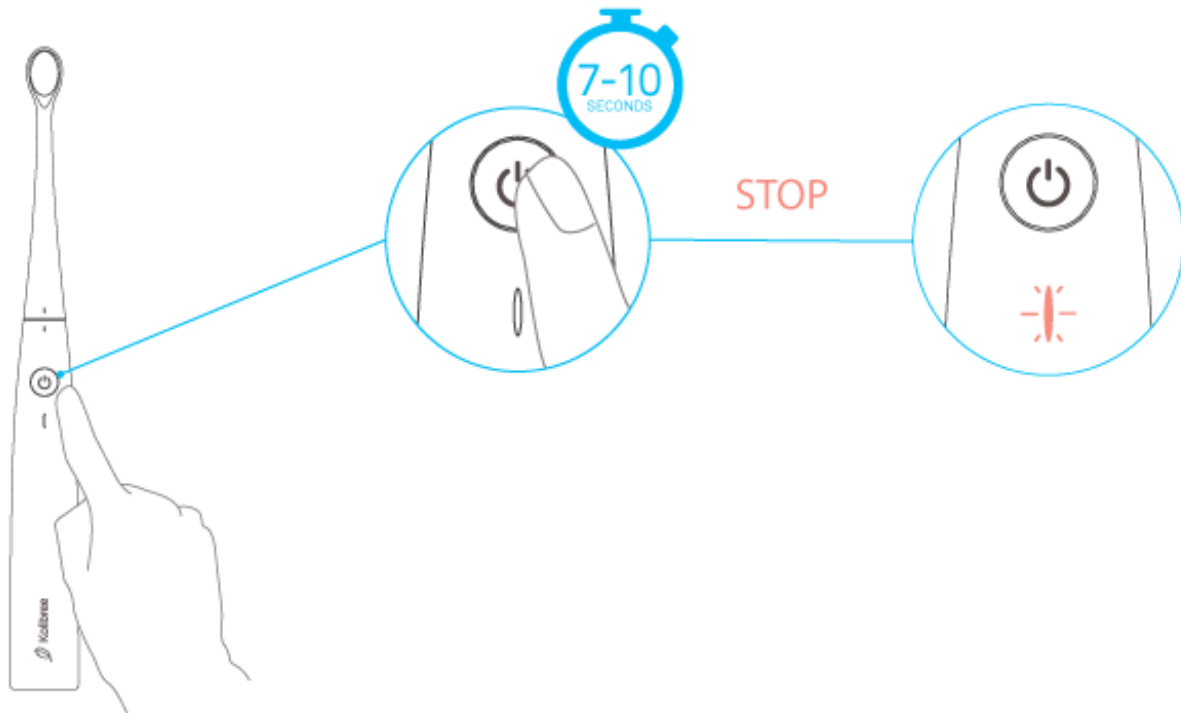


You don't have to turn your toothbrush off after every time you brush, it will automatically go into standby mode after a few minutes of inactivity.

However, if you want to turn your toothbrush off, simply hold down the button until the light blinks red. (If you aren't connected to your device then it will blink blue first and then red). At that point you can let go of the button and the toothbrush will be turned off.



Indicator lights



The indicator light located under the button on the toothbrush handle indicates different status depending on the colors displayed:

-  Green:the battery level is full.
-  Yellow:the battery level is about half full.
-  Red:the battery level is low.
-  Flashing red:the battery level is very low. The toothbrush should be charged as soon as possible
-  Blue:the toothbrush is in Bluetooth® connection-mode

2. Connecting the Kolibree toothbrush

Compatibility

The Kolibree toothbrush and Kolibree mobile app are compatible with all modern smartphones and tablets equipped with Android 4.1+ or iOS8+ systems.

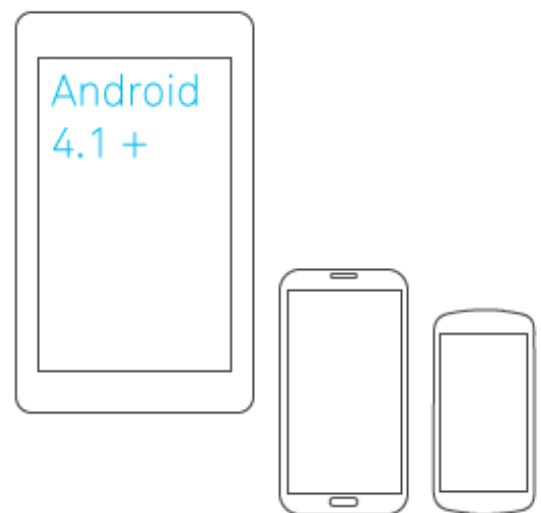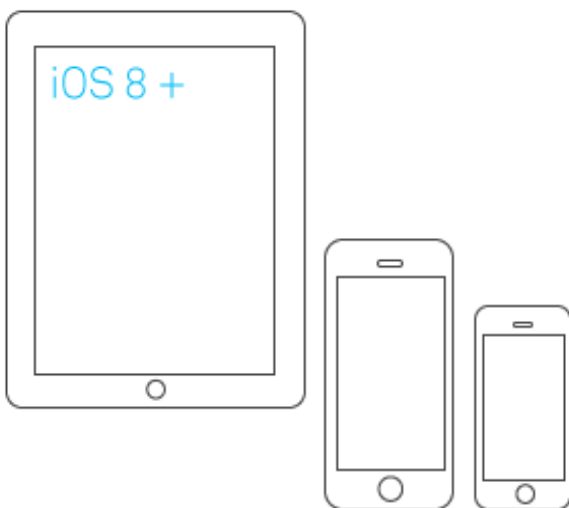


Detecting your Kolibree toothbrush with Bluetooth®

In order to connect your toothbrush to your smartphone, both the Kolibree toothbrush and your smartphone or tablet need to be turned on.

To turn on the toothbrush, press the button. The indicator light will briefly flash and then light up green.



In order to connect the toothbrush to Bluetooth®, hold down the button until the light flashes blue (pressing it for 3 to 6 seconds).

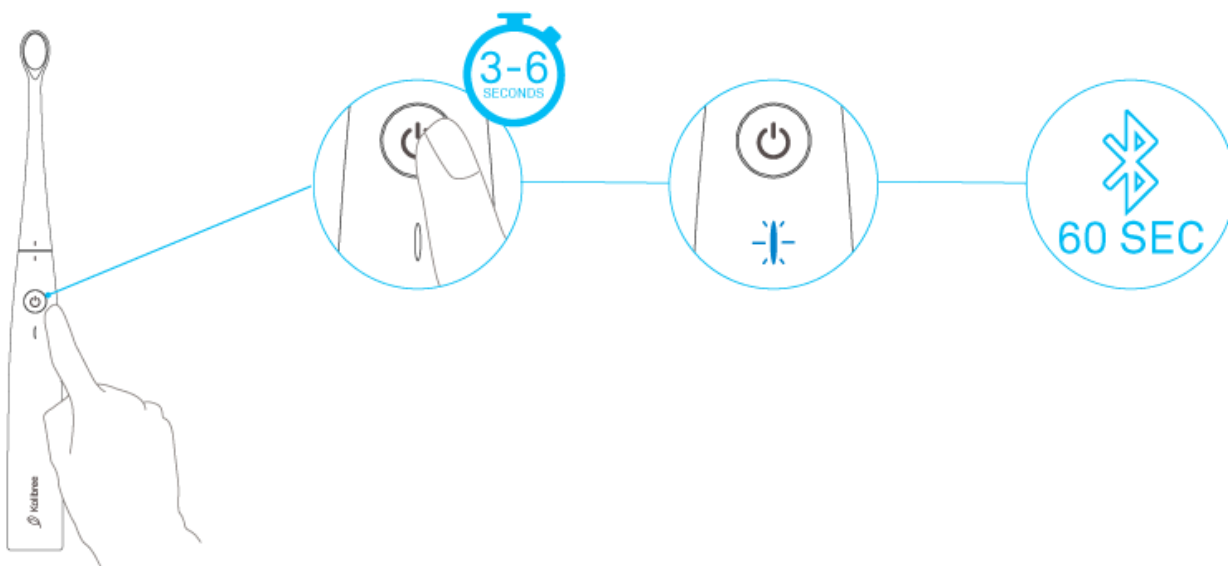

Connecting the Kolibree toothbrush to a mobile device


On your iOS smartphone or tablet, go to the “Settings” menu for Bluetooth®, enable Bluetooth® if it isn’t already enabled, and then connect the toothbrush that appears with the name KLTB001-XXXXXX, or with the MAC address that can be found on the label on the back of the box.

If you are connecting your toothbrush for the first time, you will be asked for a pairing code. Use “0000”.

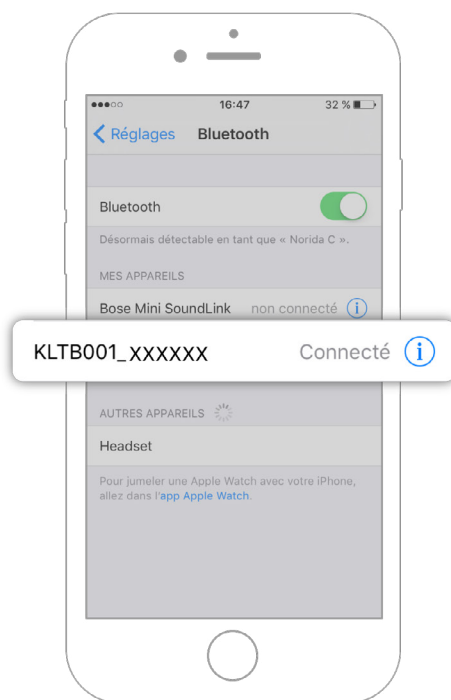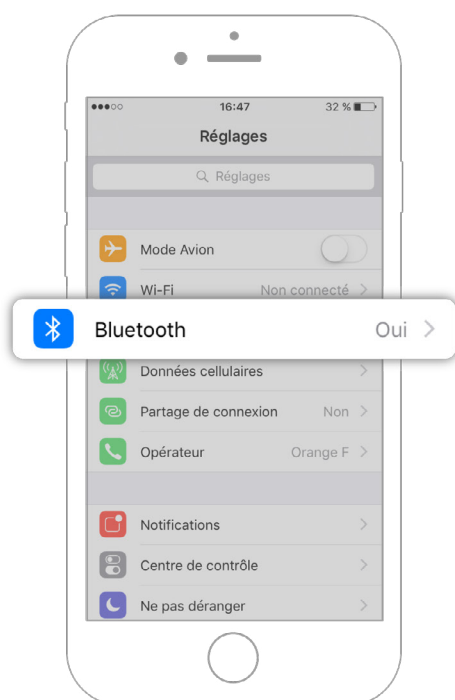


On your Android smartphone or tablet, you have two options:

- go to the “Settings” menu for Bluetooth®, enable Bluetooth® if it isn’t already enabled, and then connect the Kolibree toothbrush that appears with the name KLTB001-XXXXXX or with the MAC address that can be found on the label on the back of the box. If you are connecting your brush for the first time, you will be asked for a pairing code. Use “0000.”

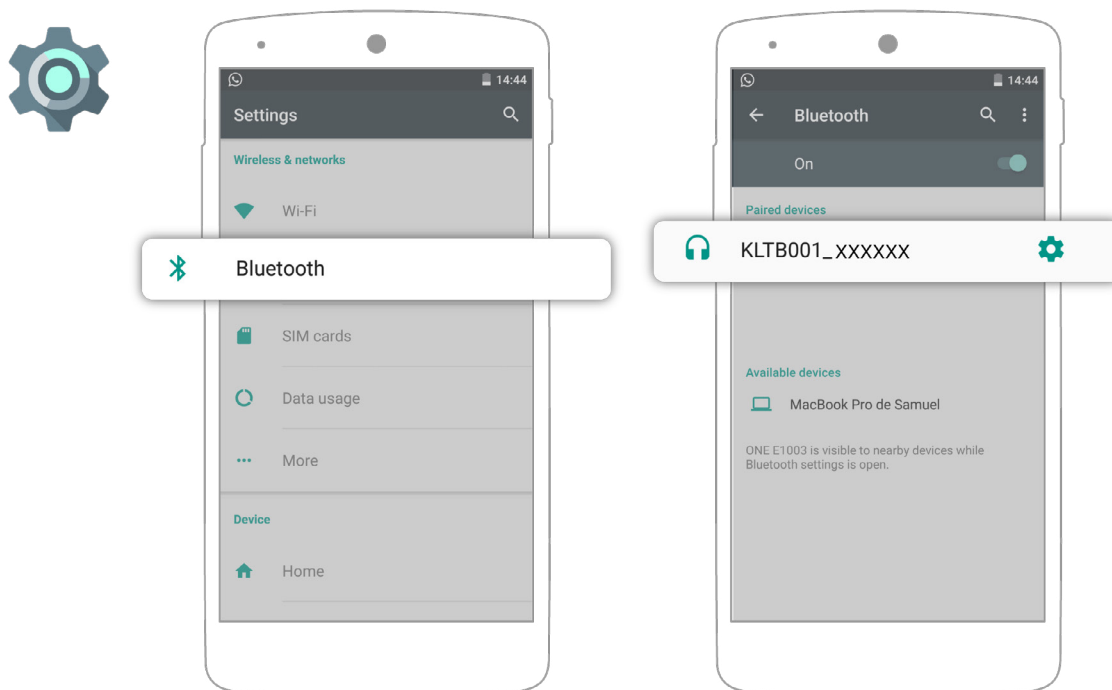

- in the Kolibree app, you can press the button showing a toothbrush at the top right with a gray cross. In the window that opens, select the Kolibree toothbrush that appears with the name KLTB001-XXXXXX or with the MAC address that can be found on the label on the back of the box.

If you are connecting your brush for the first time, you will be asked for a pairing code. Use “0000.”

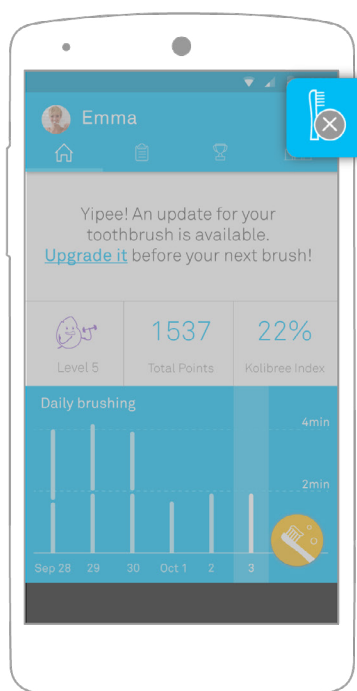

Linking toothbrushes to an account


The Kolibree toothbrush can be set to be used by a single user or by several family members. This setting is determined when the master account is created when the brush is first connected to the app.

If you choose the individual use mode, this means that the toothbrush will be used by a single person. In this case, you don't have to use the toothbrush with the app to save brushings. You can brush several times, and the data will be stored in the toothbrush until you connect to the Kolibree app again. Then the stored data will be transferred at that time.



If you choose the shared user mode, this means that the toothbrush will be used by several family members. In this case, each user must connect to his or her profile in the Kolibree app and brush with one of the activities offered. If brushing is done without the app, data will not be stored and cannot be assigned to a user.

iPhone

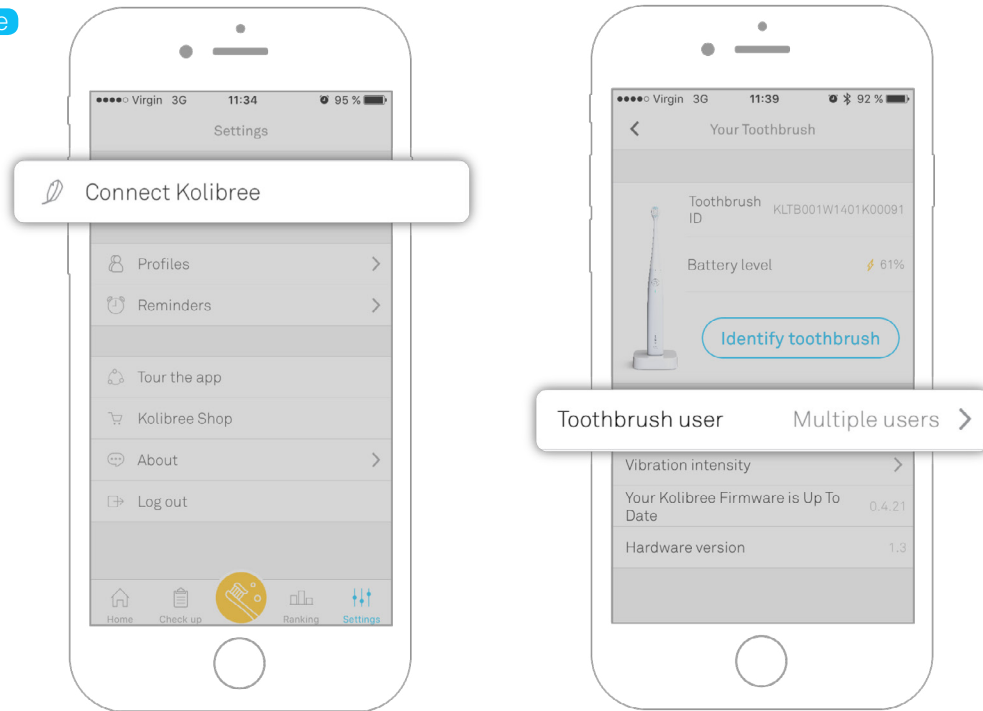

It is possible to change this setting in the app settings: Settings > Toothbrush

From the “Toothbrush” menu, you can see the connected toothbrush and define its settings.

You can change the toothbrush use mode and set it as shared or single user.

If you have several toothbrushes, you will need to connect them one at a time to change each one’s settings.

3. Brushing teeth


Changing the brush heads

Kolibree brush heads are easily interchangeable so each member of the family can use the same toothbrush, by simply changing the brush head at the appropriate time.

To change brush heads:


1. Unlock the brush head by turning it counterclockwise until the two indicators are lined up.
2. Pull the brush head upward to remove it from the handle.
3. Replace with a new brush head by lining up the indicators.
4. Turn the brush clockwise until it locks.

You can order new brush heads at any time on the [Kolibree website](#).

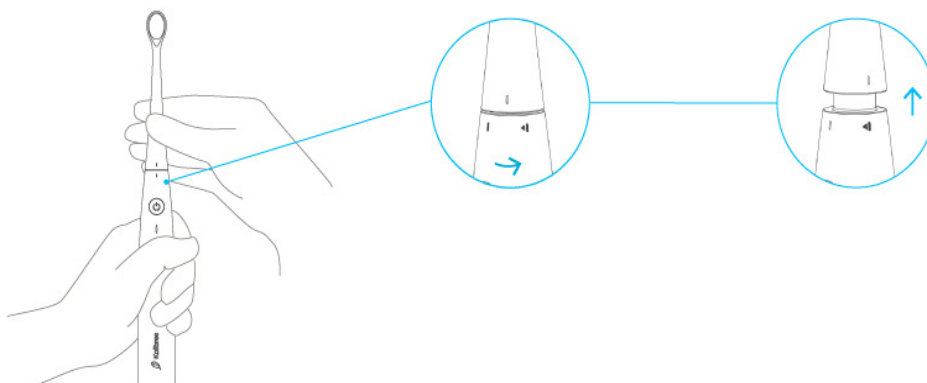


Note:

It is advisable to change the brush head every 3 months.
Use only Kolibree brush heads with your Kolibree toothbrush.

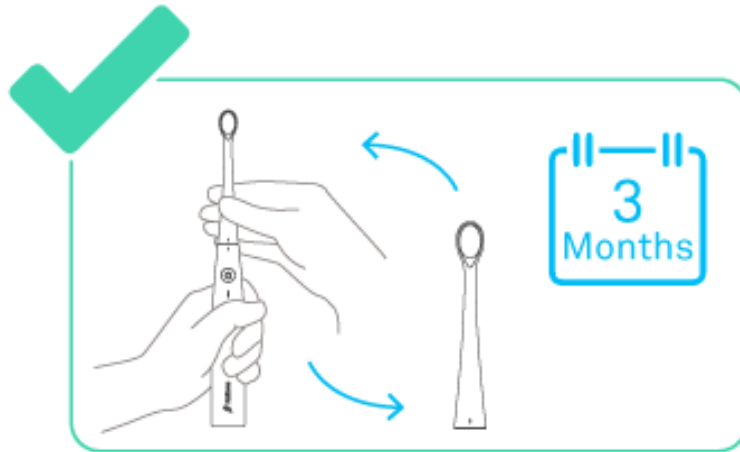

Brushing technique

It is recommended to brush teeth twice a day, in the morning and evening, for 2 minutes each time.
If you follow this recommendation, you will ensure optimal dental health.

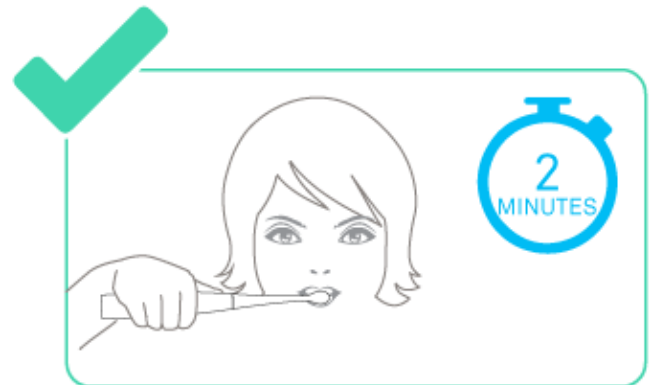

Wet the brush bristles and put a dab of toothpaste on the brush.

To avoid toothpaste being splashed or projected, place the brush on your teeth before you turn on the toothbrush.

Place the brush on your teeth at a 45° angle. Turn on the brush by pressing the button. While the vibration deep cleans, make small movements from the gums to the tips of the teeth. Clean all the outer parts of the teeth, and then the inner parts. Finally, clean the tops of the molars with small circular movements. When you are finished brushing teeth, press the button to stop vibration.

Do not press too hard on your teeth and gums during brushing.

Brushing teeth with the app


To brush your teeth with the Kolibree app activities: Turn on your toothbrush, open the Kolibree app, and select your profile.

If the toothbrush is not connected, the toothbrush icon in the upper will be gray.

-On Android: click the icon and choose your Kolibree toothbrush from the list.

-On iOS: exit the app, go to "Settings" and follow the directions for connecting.

In the list of activities, choose the one you want to use for your brushing and follow the game's instructions.

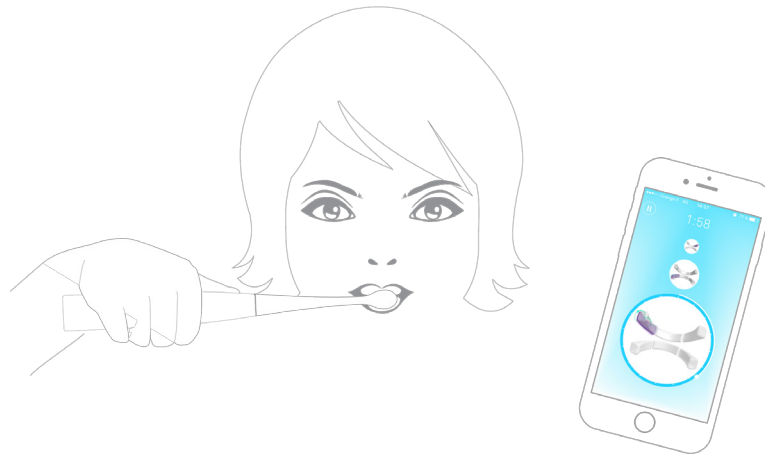

Brushing teeth without the app

If you want to brush your teeth using your Kolibree toothbrush without using the mobile app, just press the button to turn the toothbrush on or off, and use it like a regular electric toothbrush.

However, if you want to retrieve information from a brushing that was done without being connected to the app, the toothbrush automatically saves it and the data will be transferred to the app the next time you connect to it.

Note that data saving and recovery afterwards is only possible if the toothbrush has been set to individual operating mode.




4. Kolibree account and users


Installing the app on a mobile device

On your iOS smartphone or tablet, go to the App Store to download the Kolibree app and install it.

iPhone

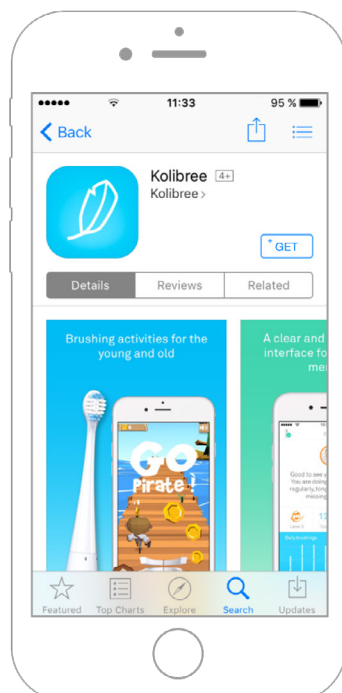


App store

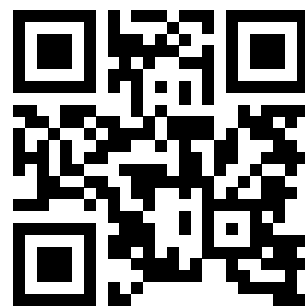


You can access it by scanning this QR code with your smartphone.

On your Android smartphone or tablet, go to the Google Play Store to download the Kolibree app and install it.

Android

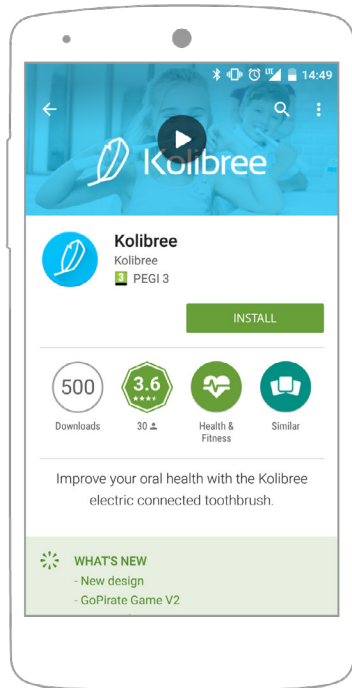

Google Play

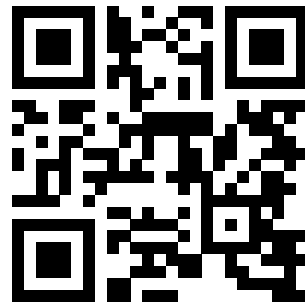


You can access it by scanning this QR code with your smartphone.

Creating the master account


If you are connecting for the first time, you will automatically be asked to create a user profile. This is the master profile that will manage your other family members' profiles if you add them.

Only the user associated with the master account can change the password and manage the other profiles linked to this account.

iPhone

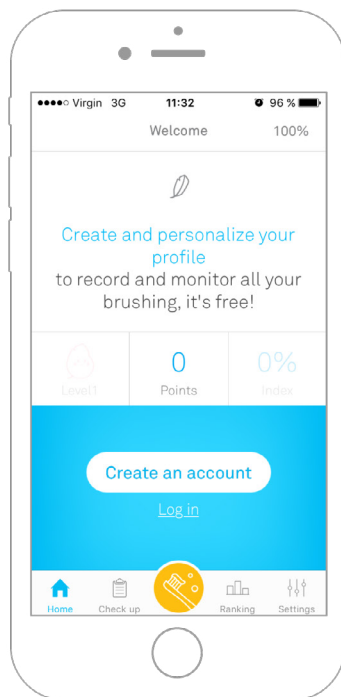

Android

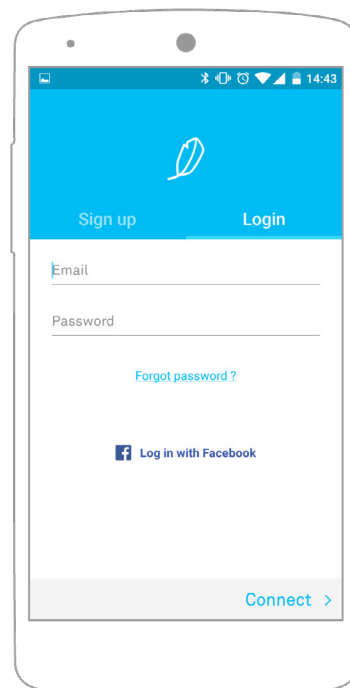

To create a new account, select “Register” on the app’s home screen and follow the instructions.

You can also create your account by linking it to your Facebook account. To do this, select “Connect with Facebook”.


You can change your personal information at any time by going to the “Settings” tab in the menu.

Connecting


If you have already created your account for the Kolibree app, you can connect in two ways:

- If you created your account with your email address, enter your e-mail address and your password.
- If you created your Kolibree account with your Facebook account, select “Connect with Facebook”.

iPhone

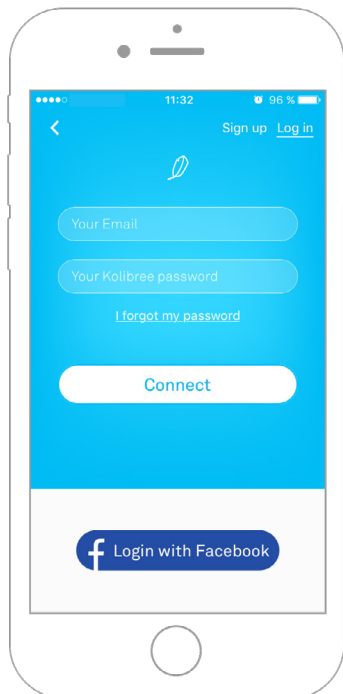

Android

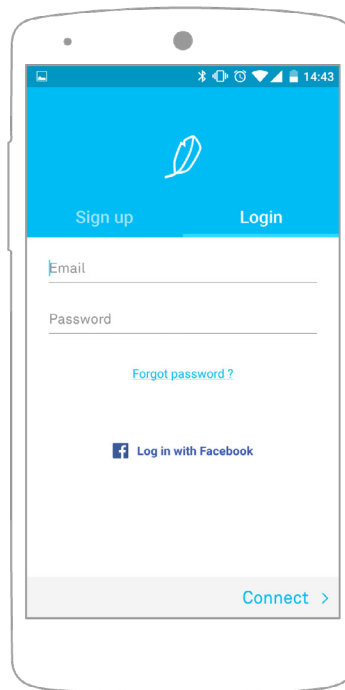


Adding users


After the master account has been created, it is possible to add user profiles, for example for each family member.

To do this, in the list of profiles, add a profile, and then follow the steps to create a new profile.

iPhone

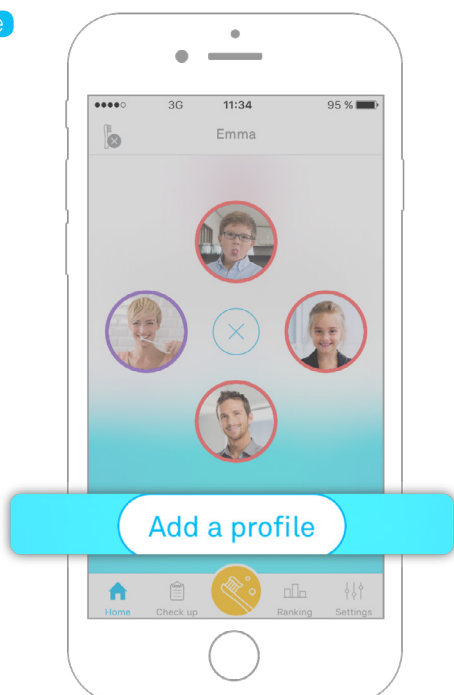

Android

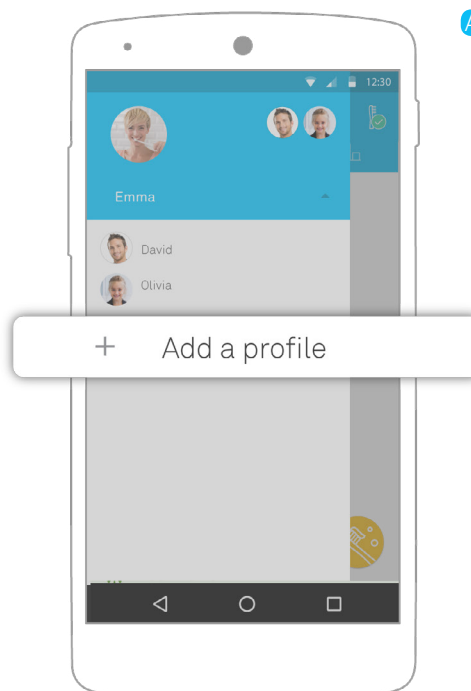


Modifying account or user information


Each Kolibree user can change the information in their profile.

To change profile information go to the “Settings” menu in the application and then choose “Profiles”. Choose the profile photo at the top of the screen by swiping left or the right until you get to the profile that you want to change.

iPhone

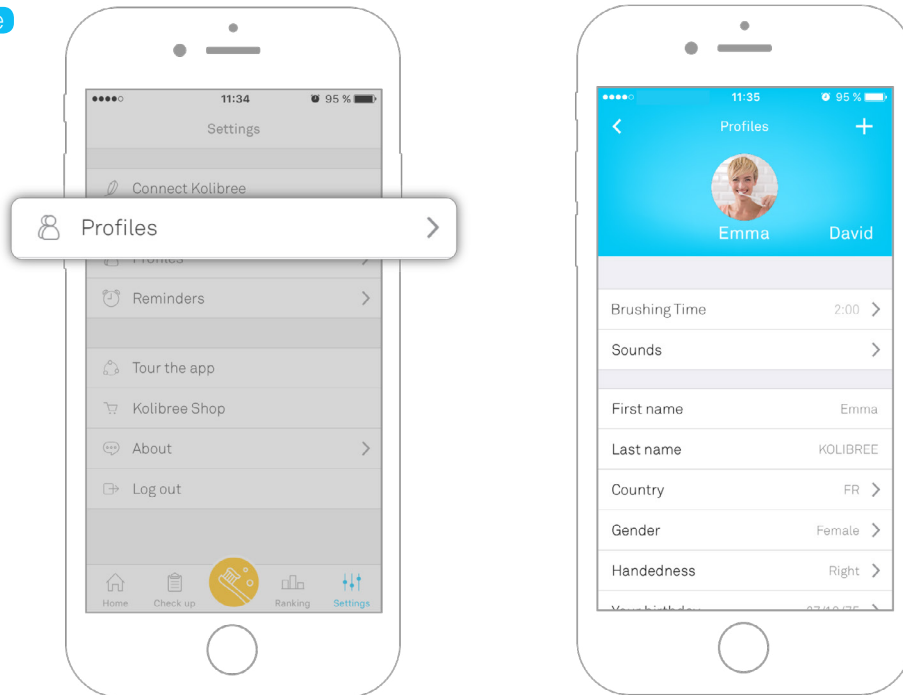


Android

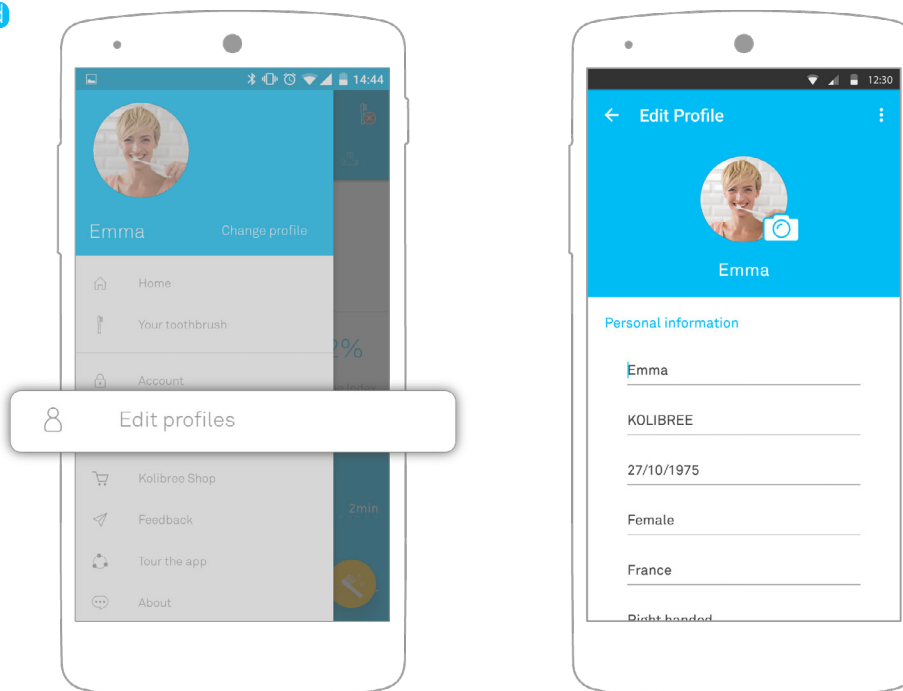


Deleting a profile

Go to “Settings” and then “Profiles”. Choose the profile photo at the top of the screen by swiping left or the right until until you get to the profile you want to delete. At the bottom, select “Delete profile”.

iPhone

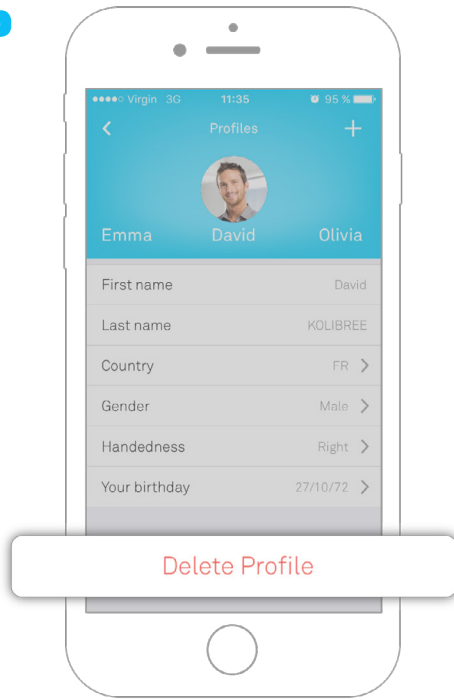

Android

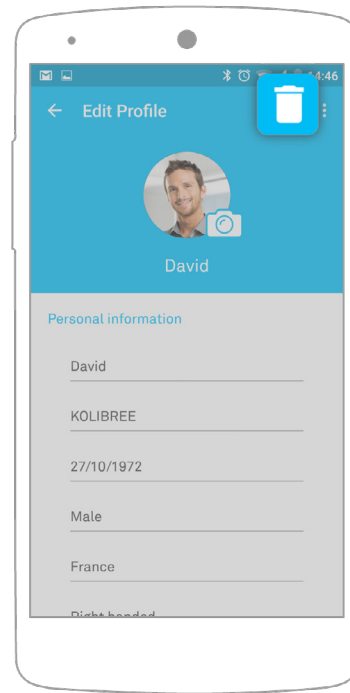

Deleting the Kolibree account

If you delete your Kolibree account, you will lose all saved profiles. Deleting an account is irreversible and will erase all data.

If you want to delete your account, please contact us.

Forgotten or lost password


If you have forgotten or lost your password, you can request for it to be reset.

Go to "Connect" and select "Forgot password". Then follow the instructions to have it reset.

iPhone



Android

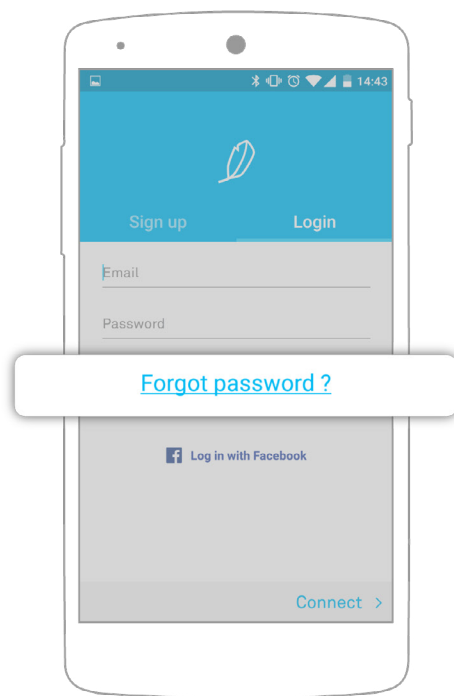


Changing the password


Only the profile that manages the master account can change the account password. Reminder: the master account is the first account created in the app.

To change the password, go to the app's "Settings" menu and then to "Account".

iPhone

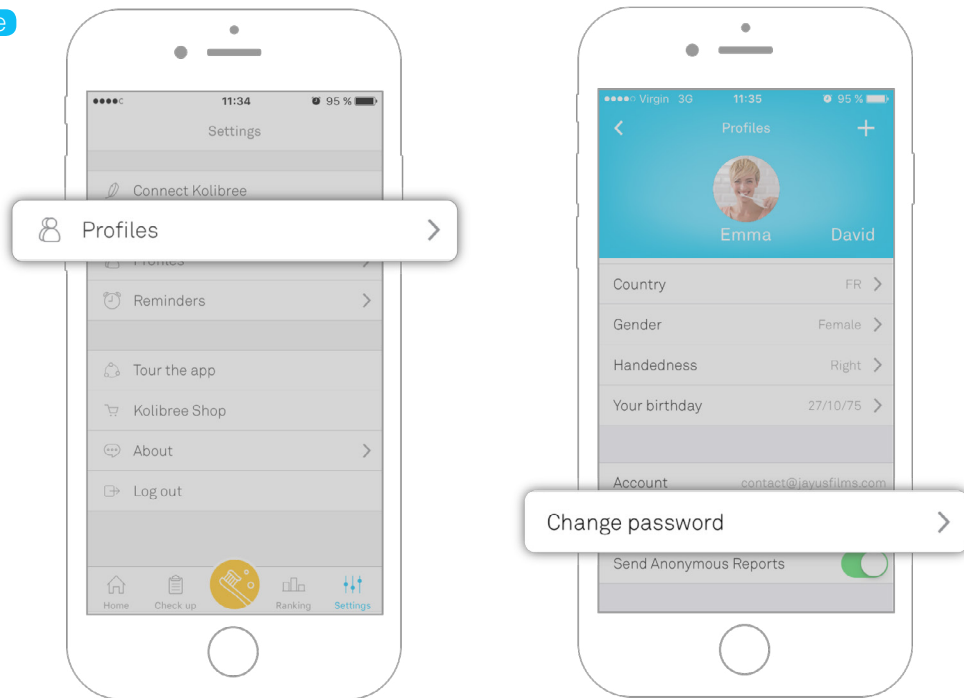

Android

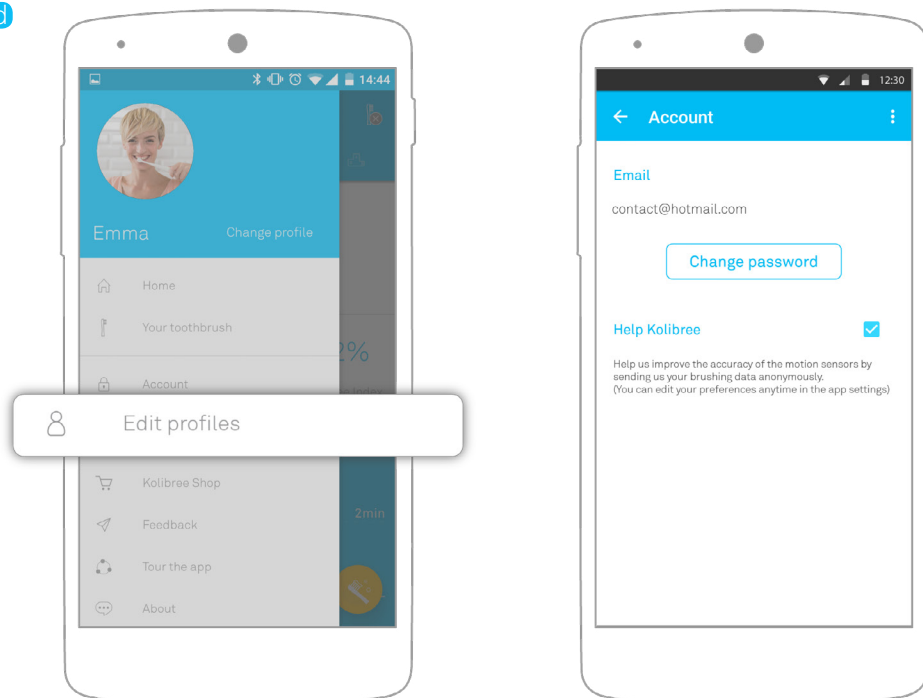


Logging out of the account

To log out of your account:


-On Android: open the menu tab in the upper left and select “Logout”.

-On iOS: open the menu to access the app settings by pressing the “+” at the bottom of the screen. In the settings, select “Logout”.

iPhone

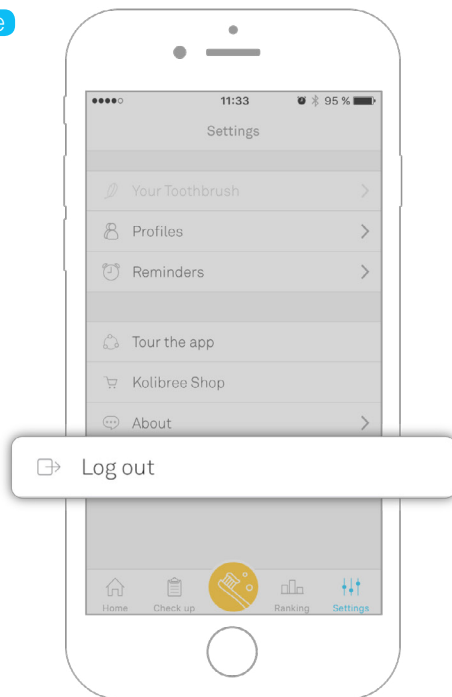

Android

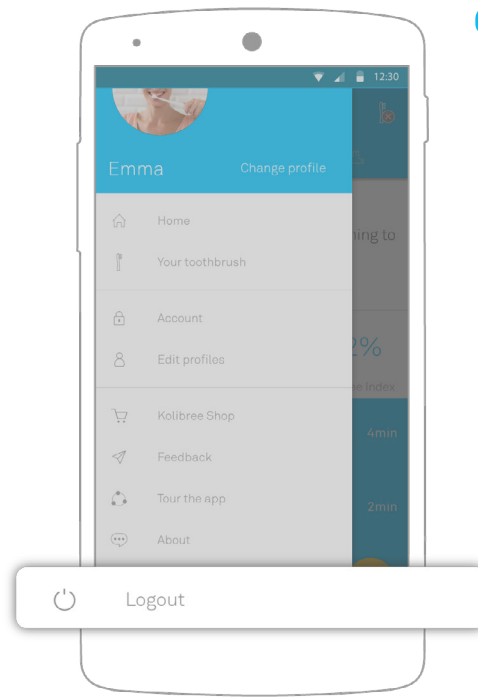


5. Kolibree App

Home


Each profile's home screen displays a variety of information

- in the center are messages from the toothbrush which vary with your activities: it could congratulate you, or encourage you to improve, and it also sends oral health tips (created by health professionals).
- your experience level
- your points
- your Kolibree index
- your last 7 days' activity

iPhone

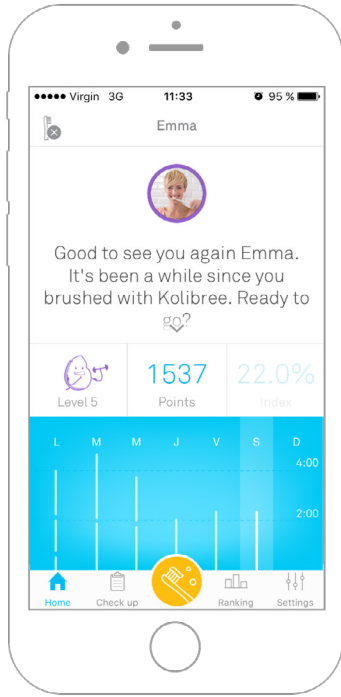



Android

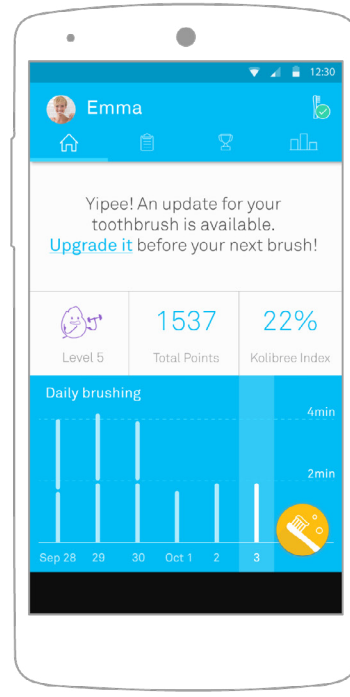


Level and badges

You can have a quick glance at your experience level in the Kolibree application. If you want more details and want to see the badges you've won, then click on your level.

iPhone

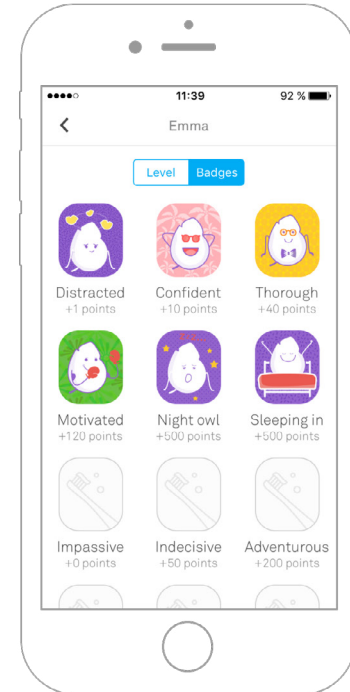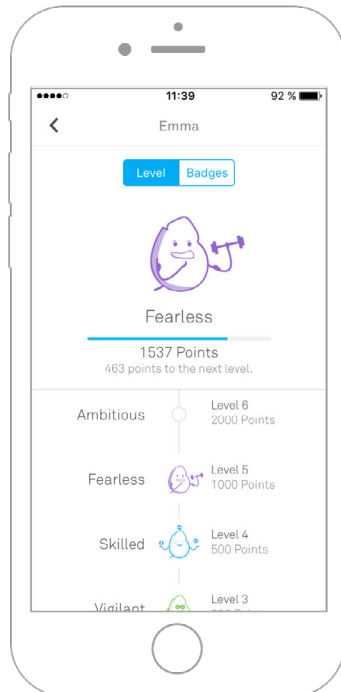

Android

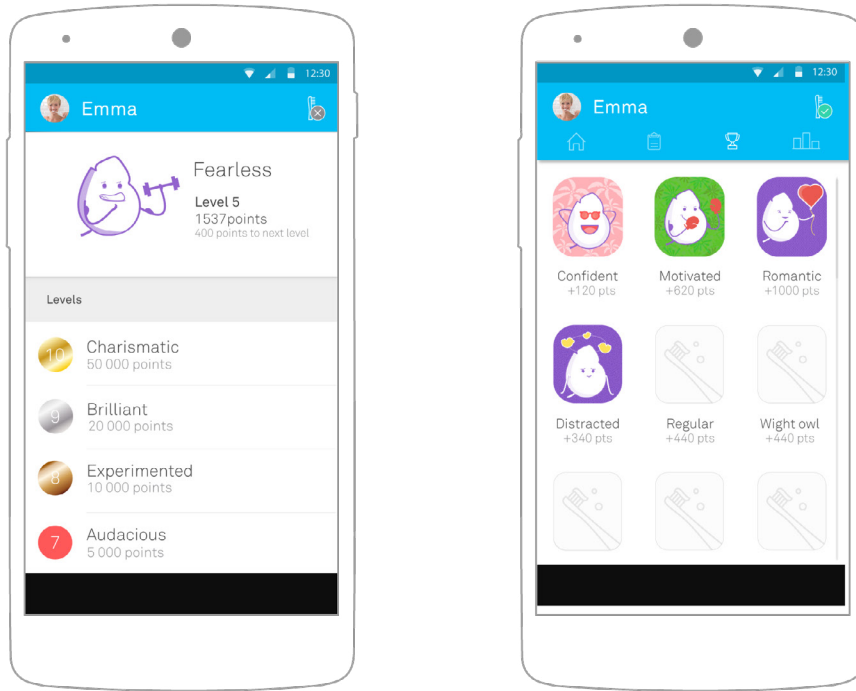


The more points you get for brushing your teeth, the higher your level. Each level has its own Kolibree cartoon mascot and achievement word. Click on the mascot to get the saying for the level. You can click on any badge you have won to get more details, or click on a badge that you don't have in your collection yet to get tips on how to collect it.

Points


Each time you brush you will earn points. You can see your points total on the main screen of the application.

If you click on your points you can access your badge collection.

iPhone

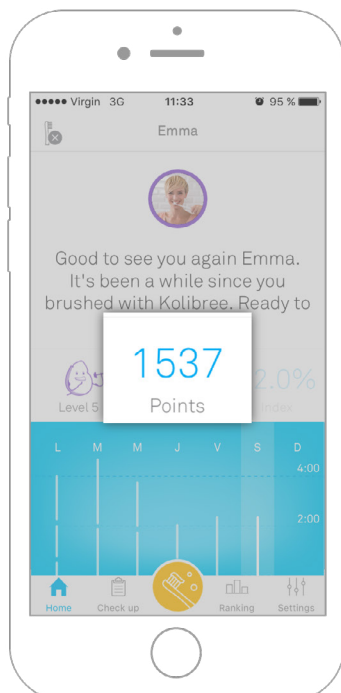

Android

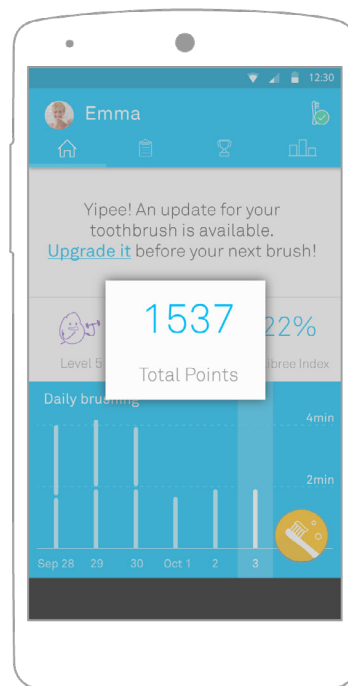

Kolibree Index

Your Kolibree Index measures the quality of your dental hygiene, according to your activity. It's based on how regularly you brush, how long you brush and the quality of your brushing.


The base Kolibree Index is 50%.

To improve it you have to


- be regular (brush your teeth every day, two times a day)
- brush your teeth for the recommended time (2 minutes each time)
- brush your teeth well (brush the correct zones for the time recommended in your activities)

Dental health professionals recommend a Kolibree Index higher than 60%.

iPhone

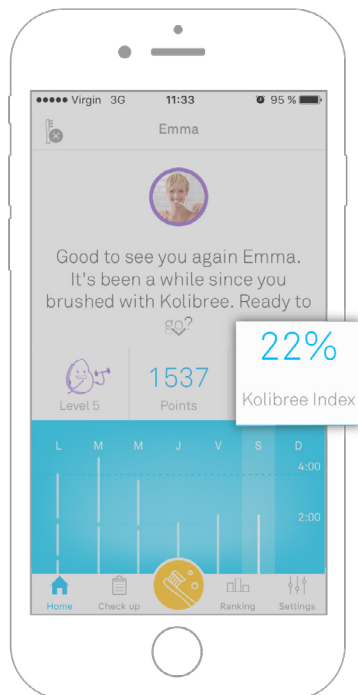


Android

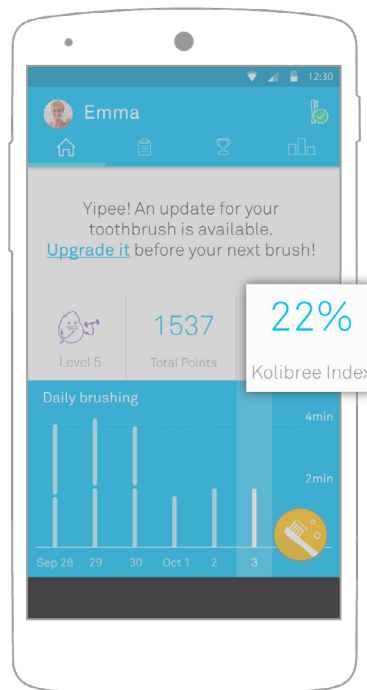


Weekly brushing

On the main screen you can check your weekly brushing on the graph. It shows the number of brushings you have done over the last 7 days. If you have done 2 brushings of 2 minutes you will see two bars on the corresponding day. Each bar can be longer or shorter, depending on how long you spent brushing.

iPhone

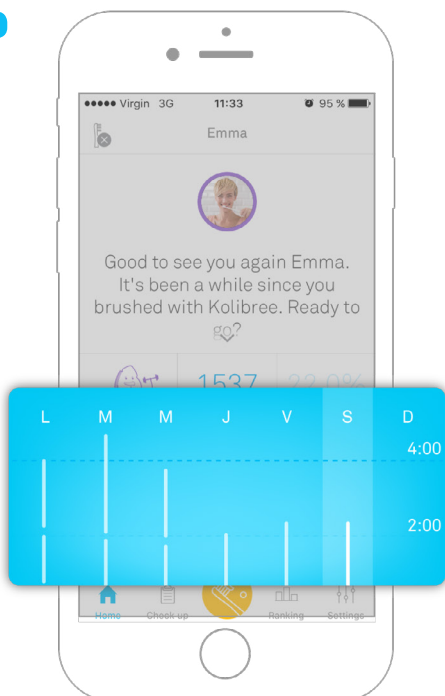


Android

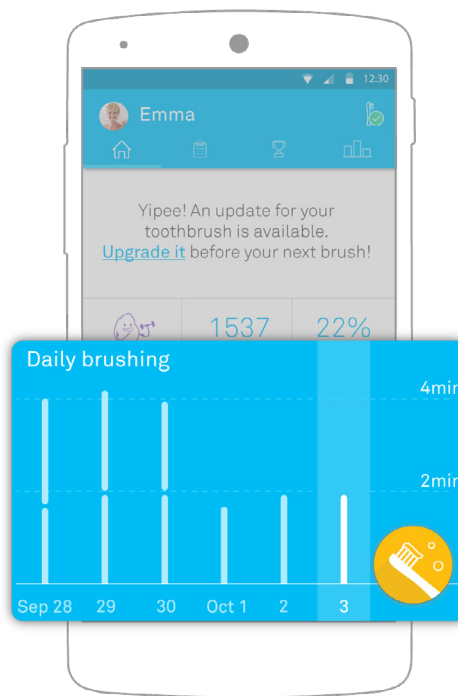


Choosing a brushing activity

Open the Kolibree application from your mobile applications. Once it's open, tap the yellow button with the toothbrush on it in the menu and choose the application you want to use while you brush your teeth.

iPhone

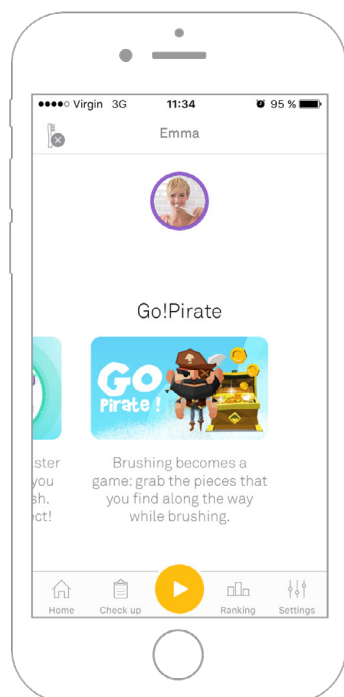


Android

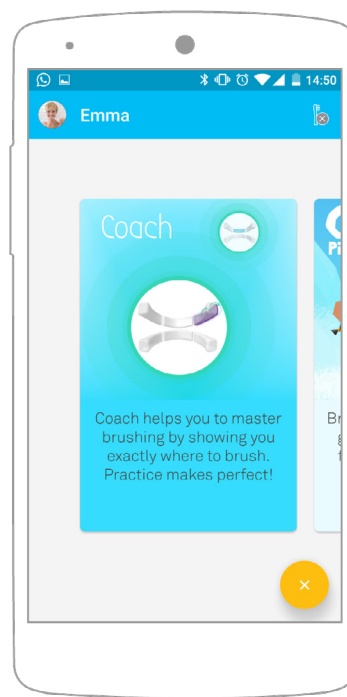


Ranking

Whatever game you use while brushing your teeth, the points you score are put in a ranking with the other members of your family who are sharing the application with you.

iPhone

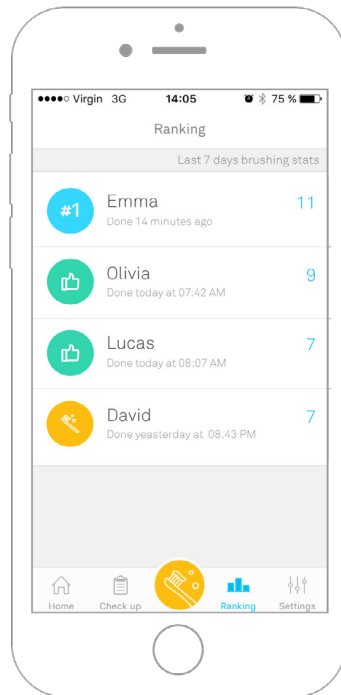

Android

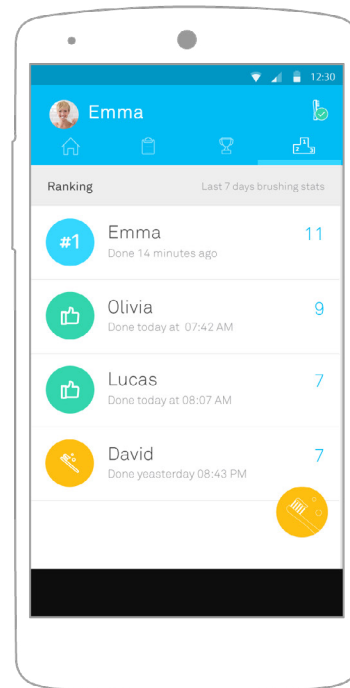

With this ranking you can keep an eye on the brushing of all the members of your family and get an overview of the last 7 days.

The icon in front of each user indicates their recent activity.


Their last brush is indicated under their first name.

The figure at the end of the line is the number of times they have brushed over the last 7 days.


Check up

To look at a more detailed version of your results go to “Check up”. Here you can see a dental diagram of the areas that you tend to neglect, based on an analysis of your last 7 days of brushing.

iPhone

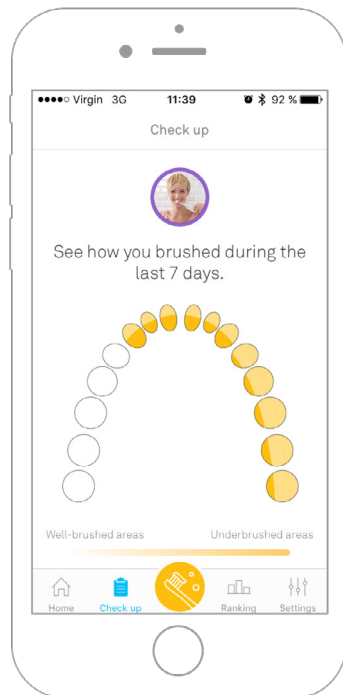


Android

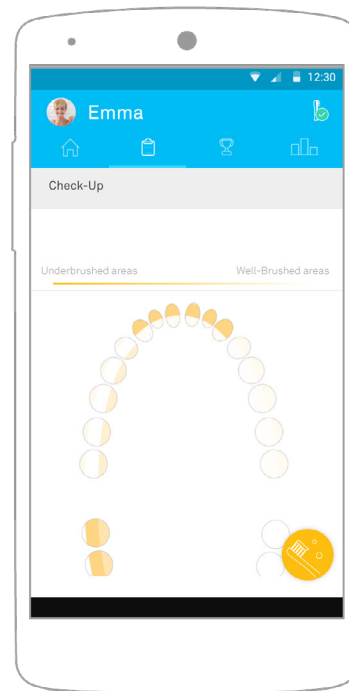


App language

The application's language depends on your phone's language. Currently, Kolibree is available in French and English, German, Spanish and Italian.

iPhone

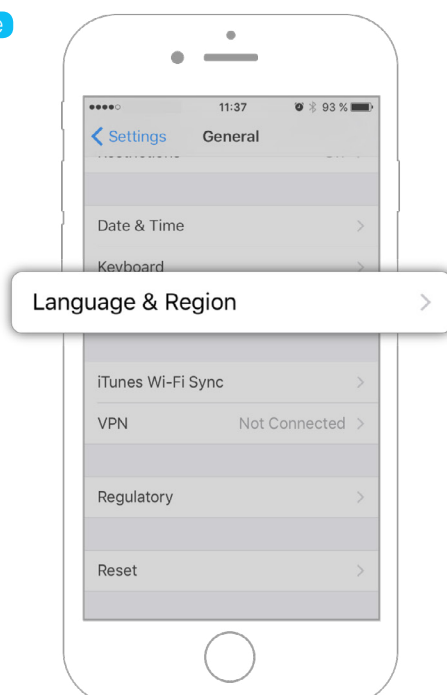


Android

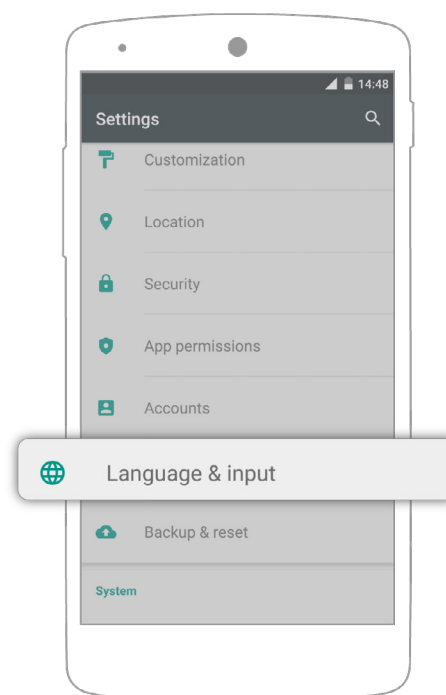


Kolibree app updates

Updates of your Kolibree app are done by your mobile device via the App Store or Google Play Store. Your mobile device will automatically offer any updates.

iPhone

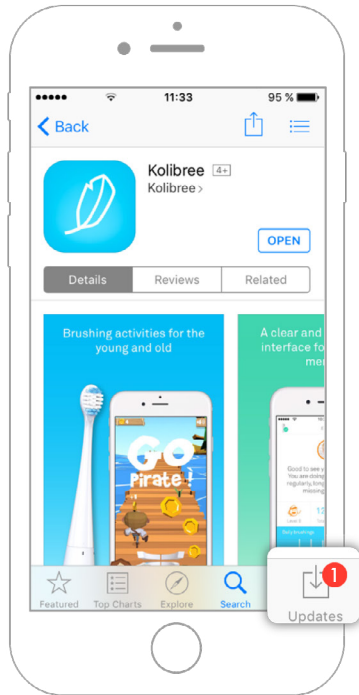

Android

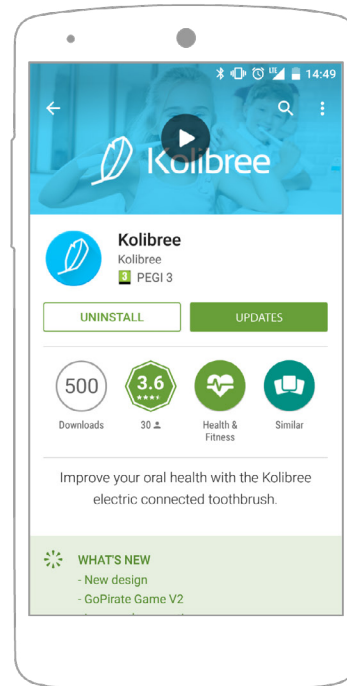


Adjusting brushing times

The default brushing time is 2 minutes.


However, your dentist may recommend that your dental hygiene requires a longer brushing time.

The default brushing time (as well as that in the Coach activity and Go Pirate! Game) can be regulated in the Kolibree app. Go to Settings > Profile > Brushing time.

iPhone

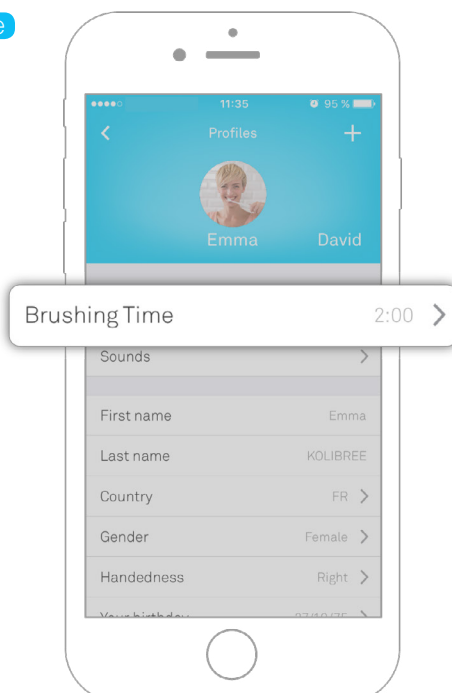

Android

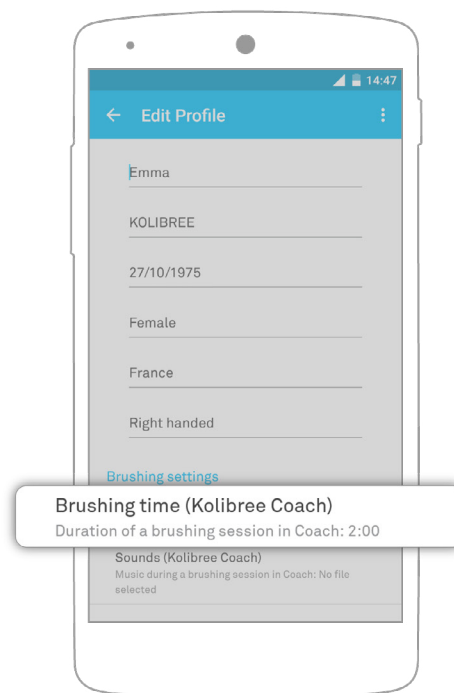

6. Brushing data

Transferring data

The toothbrush interacts with the mobile app through a Bluetooth® connection. In the same way, brushing data is transferred from the toothbrush to your mobile device.

7. Go! Pirate Game (iOS and Android)



Main

In Go Pirate! you are a pirate who needs to find Black Beard's treasure, buried on his secret island, which you can only get to once you have won the Golden Unicorn, the magic pirate galleon. Each day, thanks to your brushing, you save up so you can afford to buy a boat which will let you sail from Laguna Island to Bourbade Island and then to Corsaro Island. Finally you will once you have won enough coins to get the Golden Unicorn galleon, Black Beard's treasure will be yours!

How to play



Each time you brush, go to the zones you are shown. Your toothbrush becomes a joystick: if you brush your teeth well according to the diagrams, your character will collect coins and jewels during their trip.

Once you have chosen Go Pirate! in the Kolibree application, start the game using the button on the toothbrush.

Press the button again to pause the game.

If you want to continue brushing your teeth after pausing the game, just push on the toothbrush button again.

When the brushing game is over, the toothbrush will stop by itself without you having to press the button.

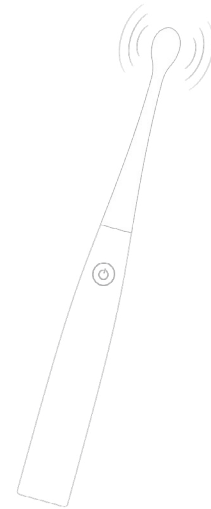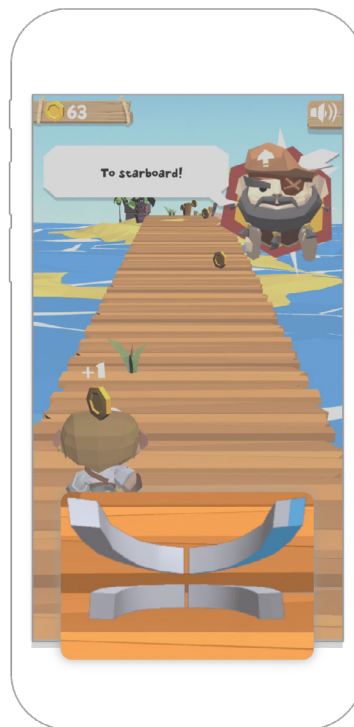

Turn on/Turn off Tutorials

Before beginning each brush, a mini tutorial will remind you to put toothpaste on the brush and then to put the toothbrush in your mouth before you push the button to start the game.


If you want to turn off these tutorials there are two ways to do it:

- when you see the tutorial before you start brushing, simply click on “Don’t show again”
- in the game settings, unclick the “Turn on Tutorials” box



Adjusting sound


You can choose to turn the Go Pirate! game sounds off. To do this, go into the game setting and choose “Sound”.

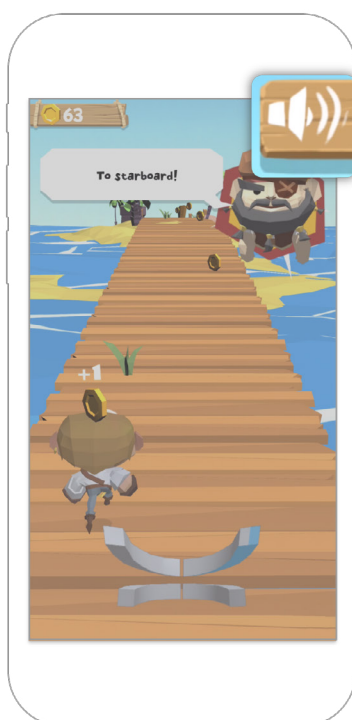


Brushing zones

During the game, Kolibree's "smart guide" will show you the areas you need to brush in real time. Brush the zones indicated by the blue light using small brushing movements over the whole surface of the teeth to really get rid of dental plaque.


In the lower levels of Go Pirate! the smart guide is divided into 8 zones:

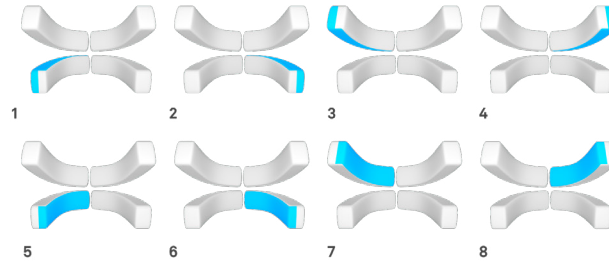

Once you get to the second world (Bourbade Island), the smart guide shows you 12 brushing zones:

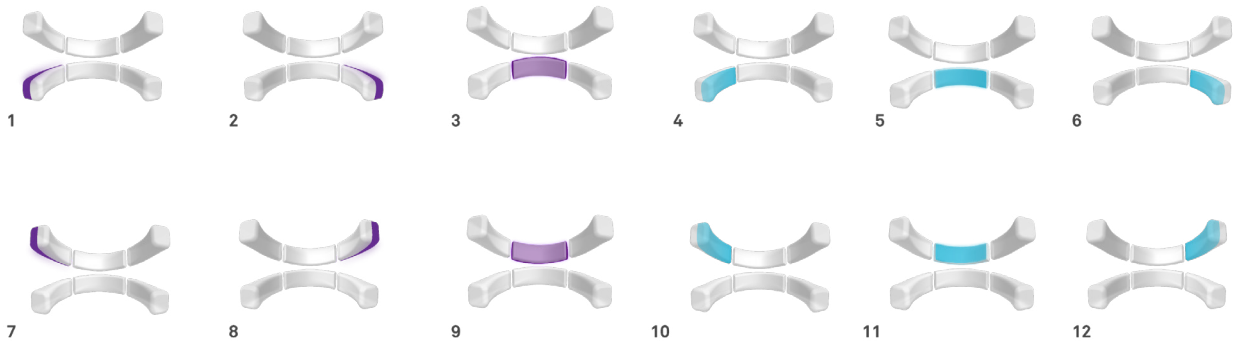


Once you get to the third world (Corsaro Island), the smart guide shows you 16 brushing zones:

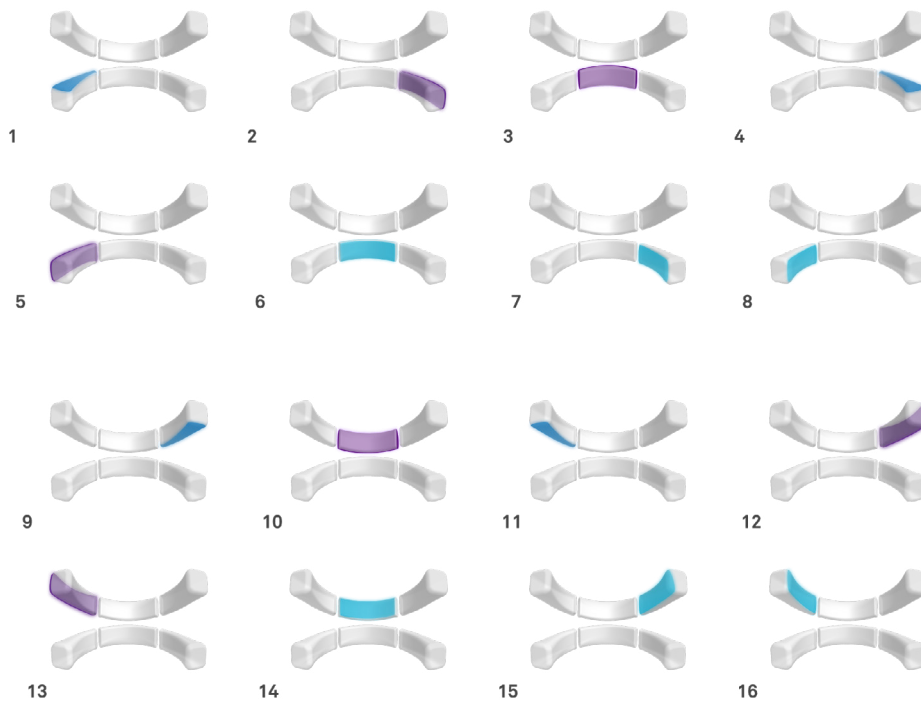

Personalize your pirate

The first time you use the app you need to personalize your pirate. If you are a girl, the application will choose a character of the same sex. Then you can personalize the skin, hair and clothing colors.

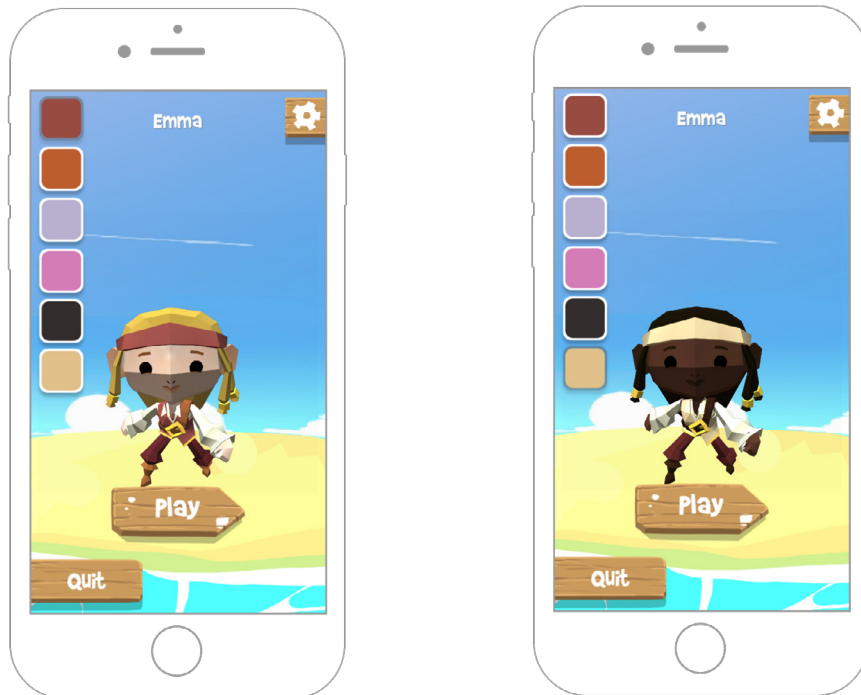

To personalize your pirate, press the setting button on the main screen and then go into the “Configure Avatar” menu.

Personalize your pirate as much as you want.

Gold coins

The objective is to get to Black Beard’s hidden island to find his treasure. During your adventure you will win gold coins every time you brush.

When you brush the zone indicated by the smart guide for the right amount of time, you will win the gold coins you find on your way.

But, once you get to level 9 (Corsaro Island) it gets more difficult: if you are not brushing correctly then the coins on your path become transparent. This means you are in the right zone in your mouth but you are not brushing the right surfaces. If you adjust your brushing so it follows the smart guide then the coins will reappear!

While you brush you can also pick up jewels that will be changed into gold coins.

The jewels are worth the following:


Emerald = 5 gold coins

Sapphire = 25 gold coins

Ruby = 50 gold coins

To get to Bourbade Island you need 3,000 gold coins.


To get to Corsaro Island you will need 24,000 gold coins. Finally, you will need 85,000 gold coins to get the Golden Unicorn galleon and go to Black Beard’s secret island.

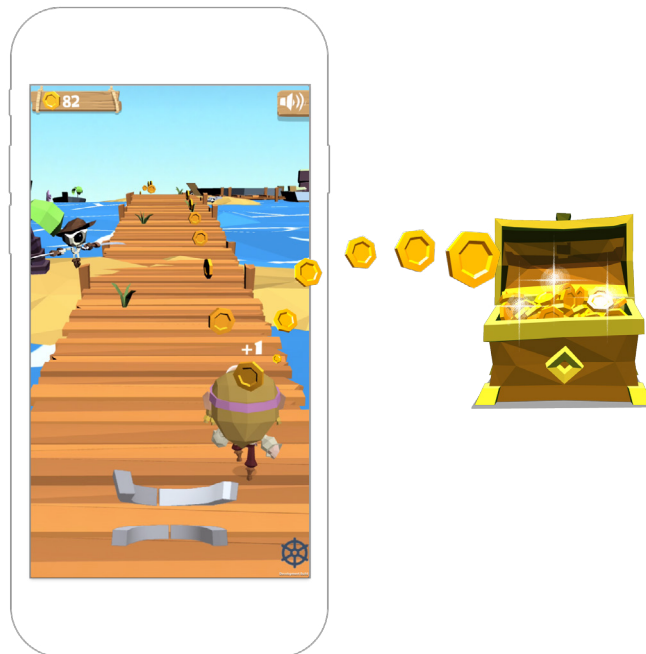

Grades

In Go Pirate! You are a pirate looking for Black Beard's treasure. You begin the adventure as a Cabin Boy but as you continue to brush you will gain experience and get promoted: sailor, helmsman, gunner, etc. Until you eventually become the Captain.

Your promotion depends on the number of brushings you do. For example, you are promoted to seaman after 10 brushings, to petty officer after 100 brushings, etc.

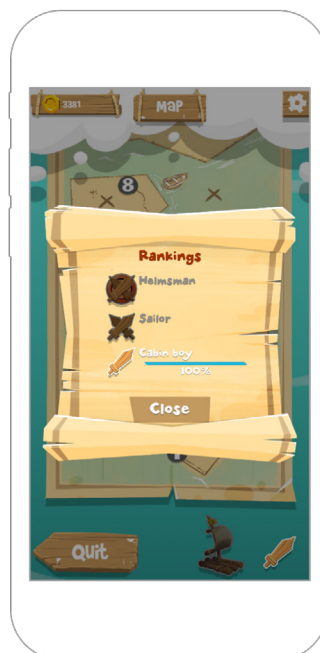

Vidéo de présentation du jeu Go Pirate!

La première fois que vous arrivez dans le jeu Go Pirate! Une vidéo de présentation joue automatiquement pour vous expliquer l'objectif du jeu.

Si vous souhaitez la revoir une nouvelle fois, depuis le bouton paramètres du jeu, allez dans "Rejouer intro".

The treasure map

On the treasure map, which you can see on the home page of the game, you can see how much further you have to go before you arrive at the secret island where Black Beard has hidden his treasure. You have to cross the seas and the islands to get there.



The steps on each island correspond to a series of brushings which will let you win gold coins so you can buy a boat to take you to the next island. You may also discover hidden treasure along the way.

Change the language

Via the game settings button you can change the language. Go Pirate! is available in English, German, Spanish and Italian.




8. Kolibree Coach Activity (iOS et Android)

Main

Kolibree Coach is the app that will help you effectively improve your oral health day after day by guiding you through a method approved by oral health professionals. Kolibree Coach shows you if your brushing is correct or not in real time.

Brush your teeth, following the areas of the mouth indicated by the app. Brush for the indicated amount of time and you will improve your dental hygiene with each brushing. The indicator will turn green if you brush correctly and red if you do not brush in the right place.

Once brushing is done, if you have brushed for 2 minutes, you'll earn coins that show your progress along with the other members of your family in the Kolibree app.

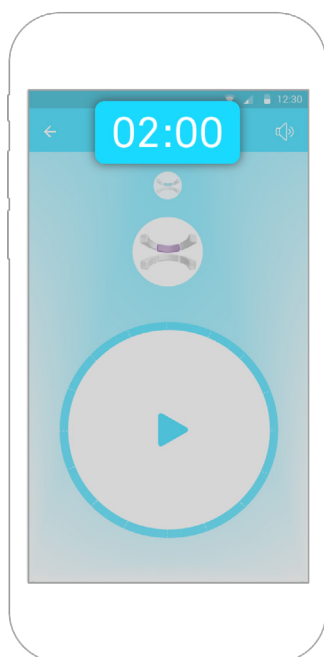


How to brush

Kolibree Coach is a tooth brushing guide which tells you how long to brush each zone so that you get to the minimum recommended brushing time of 2 minutes. Your toothbrush becomes a sort of «Dental GPS» and can detect whether you are brushing the zone that you are shown on the screen.

In this way, Kolibree Coach will help you improve your brushing technique by giving you instant feedback on how you are doing.

If you want your teeth to be well-brushed then you need to do it for a minimum of 2 minutes. Your mouth is divided into 16 brushing zones during this brushing time.


If you brush the indicated areas correctly for the right length of time then you will win points and improve your position in the Kolibree app.

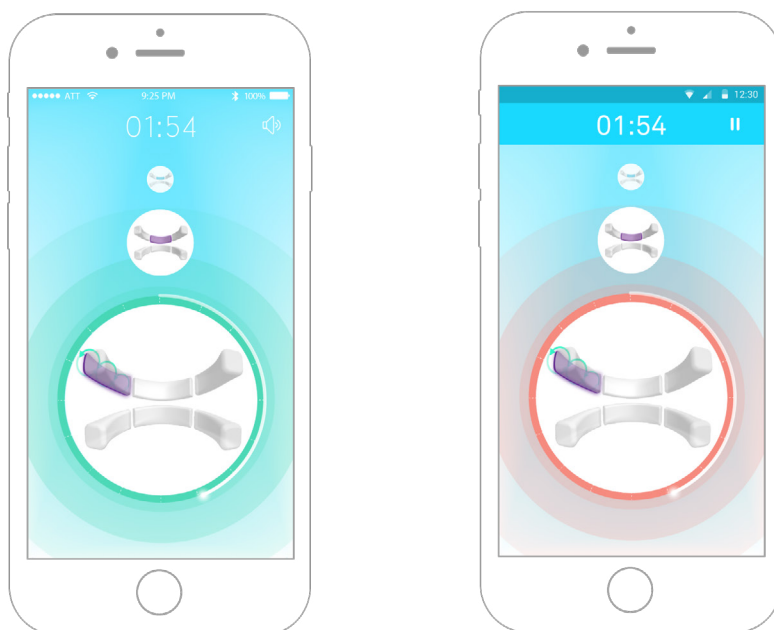

Adjusting sound


At any time during your brushing, you can turn the sound on or off by selecting the sound icon.

If you want to use your own music, you can pick a song from your mobile device's library by choosing the sound icon in the upper right.

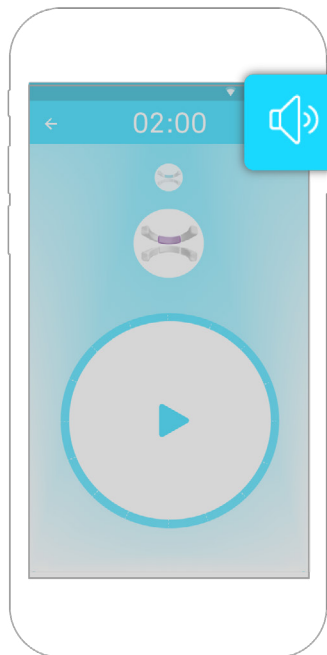


Brushing zone

The zones that are currently being brushed appear in the main center bubble, while those that need to be brushed next appear in smaller bubbles at the top of the screen. We have divided the mouth into 16 zones for the Kolibree Coach.

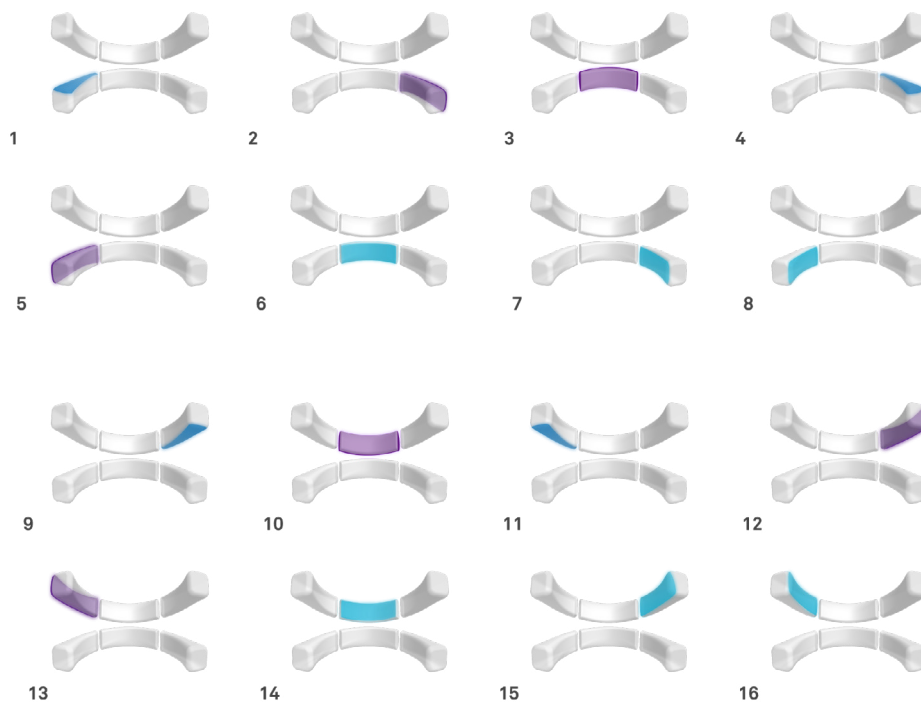


9. Kolibree toothbrush maintenance


Precautions for use

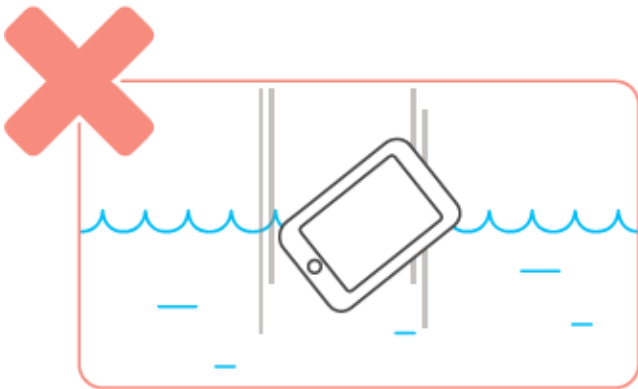


Pay attention to your mobile device during brushing to make sure it doesn't slip and drop.

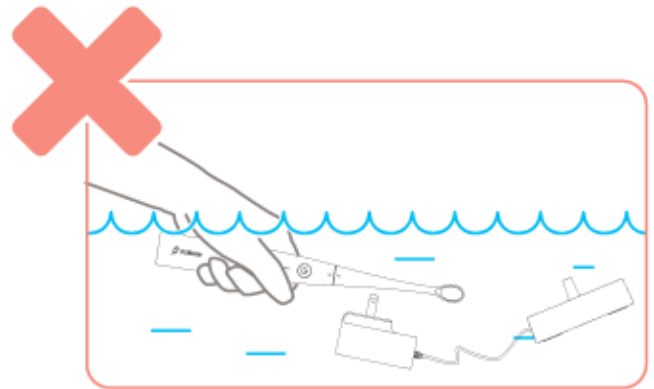


Do not immerse the brush or the base in water.

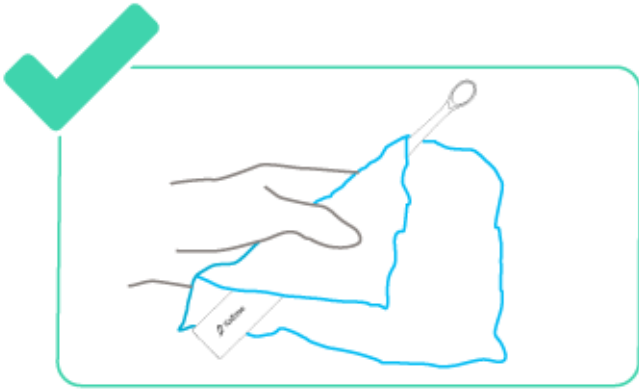


Clean the brush handle with a damp cloth.

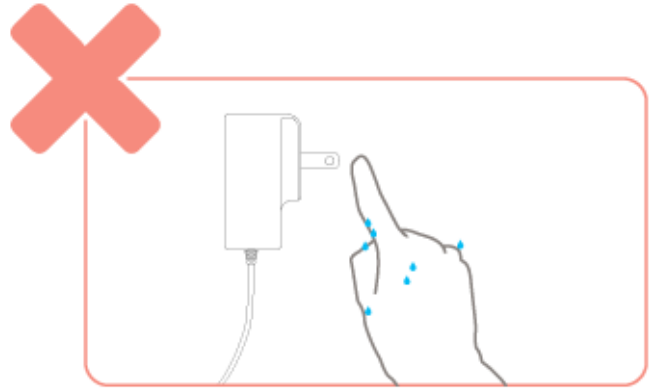

Do not touch the electrical outlet with wet hands.

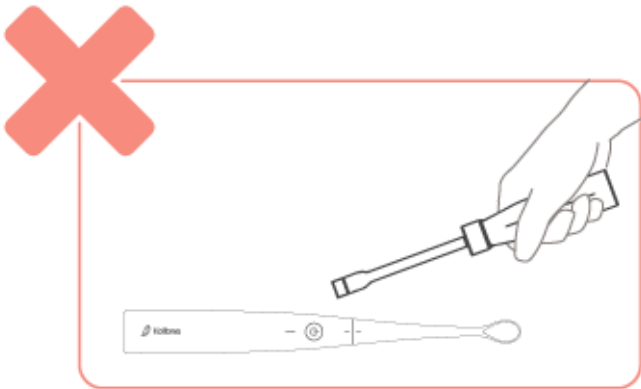





Do not modify the product.
Do not attempt to disassemble or repair.

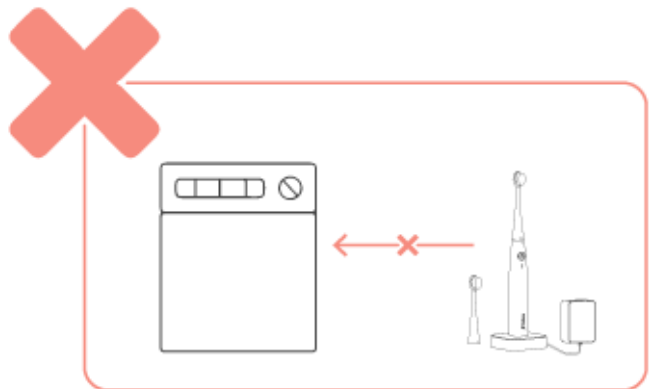

Do not clean brushes, the brush handle, or the induction charger in the dishwasher.




Cleaning brushes

After each use, clean your brush thoroughly.
To do this, you should remove the brush head from the handle and rinse it with clean water.
Then dry it before returning it to the base.

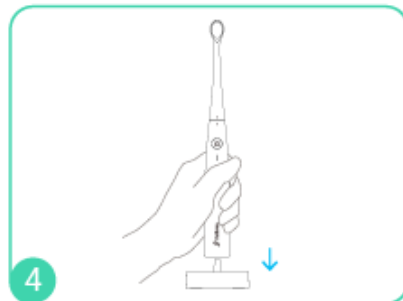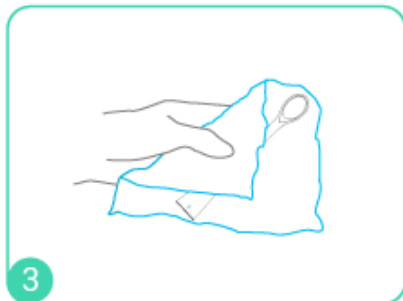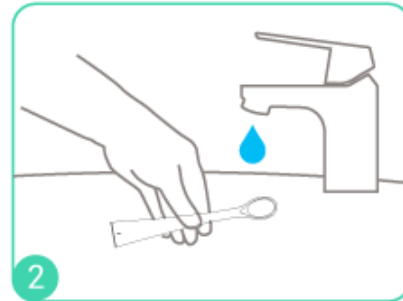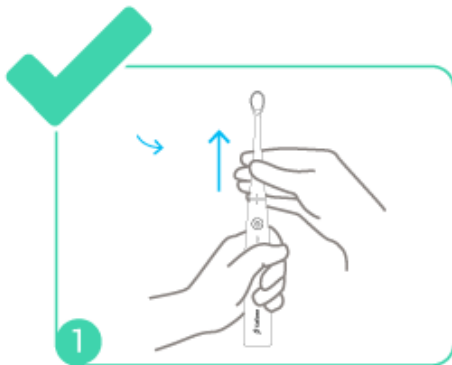

Changing brushes

Turn the brush head until the two indicators are lined up.
Pull the brush head off and replace it with a new one.
You should change the brush head every 3 months.

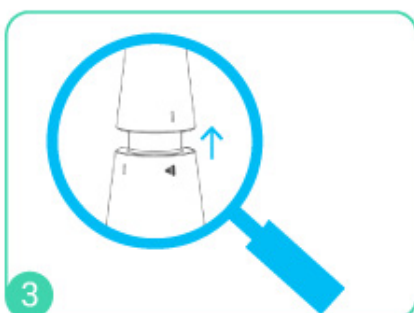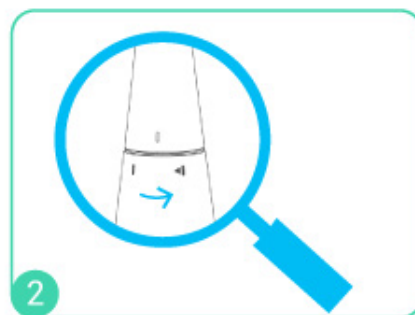
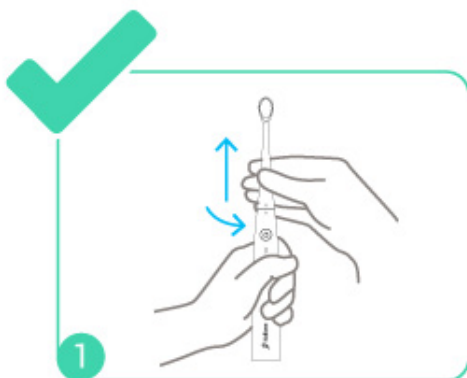


Ordering replacement brush heads

You can order your replacement brush heads on the Kolibree [online store](#) or contact the nearest authorized Kolibree retailer.

Cleaning the charger base

Regularly clean the Kolibree toothbrush handle and charger base with a damp cloth. The toothbrush handle is waterproof but rinsing it in water is not recommended. Do not clean the toothbrush or any of its accessories in the dishwasher.

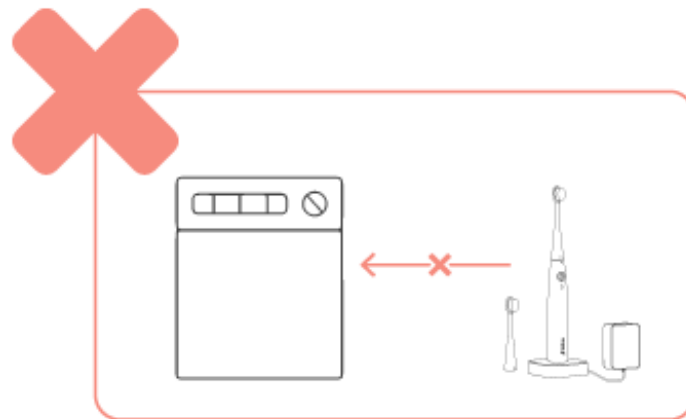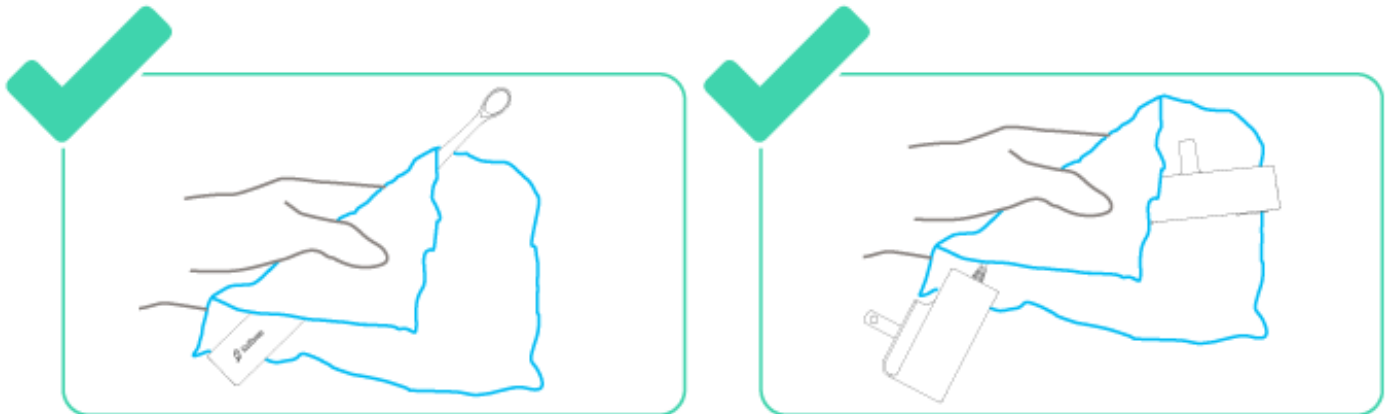





Storing the toothbrush

If you aren't planning to use the toothbrush for a long period, we recommend unplugging it from the electrical supply, cleaning the handle and brushes, and storing them in a place protected from light and heat.

Recycling the product

This product contains rechargeable batteries. To protect the environment, do not dispose of the used product with the household trash. You can place it in appropriate collection sites that are available in your country.

When the product can no longer be used, remove the toothbrush battery as follows:

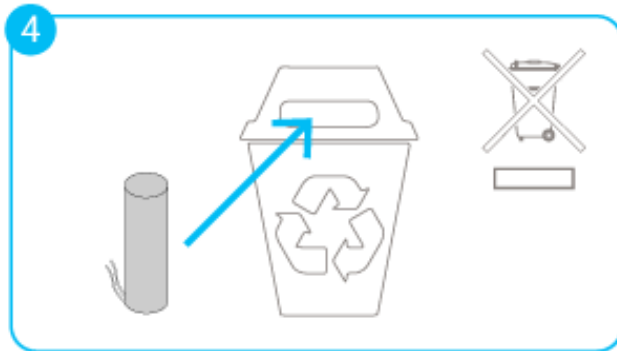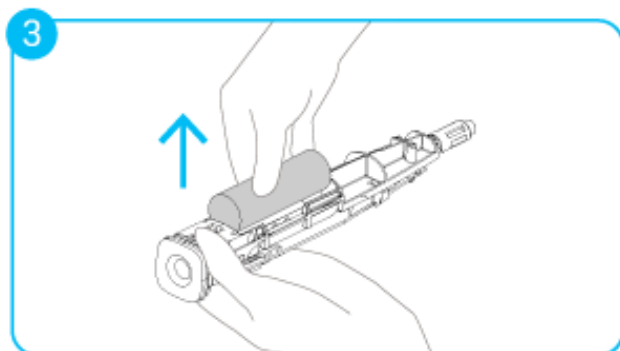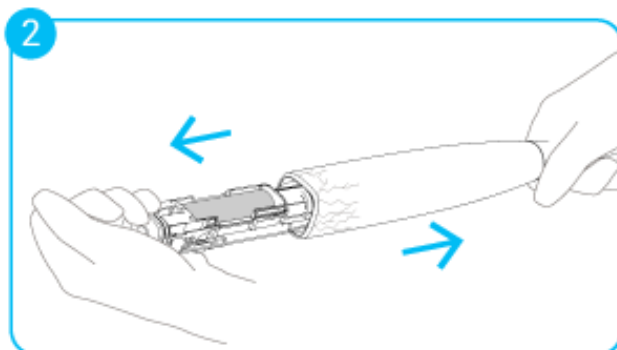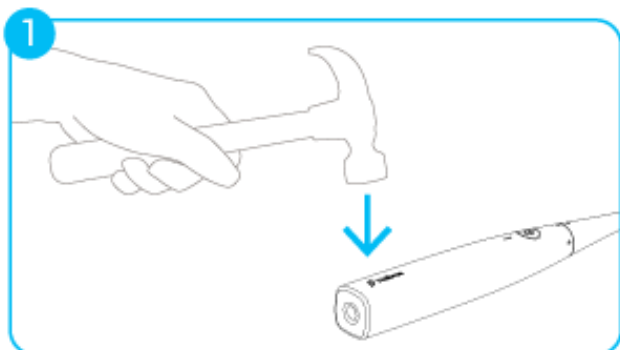

Warning! This action is irreversible. Opening the handle makes it unable to be used and invalidates the warranty. Dispose of it in accordance with the environmental regulations in force in your area.

10. Technical information

Characteristics

GENERAL	
Guarantee	1 year or the current legal period from the date and place of original retail purchase
Color	Pure white
Size	23.3 cm (9.2 inches)
Weight	70 g (2.5 oz)

TOOTHBRUSH	
Internal motion sensors	3 axis accelerometer, 3 axis gyroscope, 3 axis magnetometer
Battery	Rechargeable 650mAh Li-Ion
Autonomy	Around 1 week of typical usage*
Connection	Bluetooth® 2.1 + EDR Class 2 – 2.4 GHz
Ergonomy	Soft touch handle with a single push-button for simplicity and easy handling
OS support	iOS 8+, Android 4.1+
Vibration	15 000 vibrations/minute
Waterproofness	Water resistant
Display	one RGB LED informing on battery level and indicating the state of communication with the connected device

CHARGER	
Input	100-240V, 50/60Hz 0.3A
Output	7VDC - 500mA 3.5W LPS, induction charging system

BRUSH HEAD	
Length	9.5 cm (3.75 inches)
Bristles	Nylon bristles with rounded tips to protect gums
Material	Certified food grade plastic
Replacement pack	3 brush heads in three colors – Pack size : 15 × 8 × 1.8 cm (5.9 × 3.1 × 0.7 inches)

KOLIBREE PACK	
Content	1 Kolibree toothbrush, 1 extra brush head, 1 charging base, 1 user manual
Box size	29.1 × 10.7 × 7.8 cm (11.4 × 10.7 × 3.1 inches)
Box weight	640 g (22.5 oz)

* Two brushings of 3 minutes a day with Bluetooth® on

11. Important precautions

Dangers

Dangers

In order to reduce the risk of electrocution:

- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not use the device while taking a bath or shower.

Warnings

To reduce the risk of burns, electrocution, fire, or injury:

1. This product is not for mentally or physically disabled children or adults unless they are supervised by a responsible person. Do not let them use or play with the device alone.
2. This product is designed to be used for the purpose and in the manner described by Kolibree. Do not use any accessories that are not recommended by Kolibree.
3. Do not touch the charger when your hands are wet, or you risk electric shock.
4. Keep the charger base and electric cord away from heated surfaces.
5. When you unplug the charger, pull on the adapter and not the cord.
6. Do not do drop the device or insert any object into any of its openings.
7. Do not modify, disassemble, or repair the device under risk of fire, electric shock, or injury.
8. Do not use the device if it is defective, if it does not work properly, or if the charger is damaged or has been dropped in water. If any of these occur, contact Kolibree or your authorized Kolibree retailer to have the product examined or repaired.
9. If the toothbrush has fallen, you should change the brush head before using it again even if no damage is found.
10. The device must be operated only under very low safety voltage consistent with the unit's label.
11. This device is supplied with a charger and an electric plug included. The device must be used only with the power supply provided. Do not replace or damage any parts or you risk electric shock. Occasionally check that the charger is not damaged. Damaged or malfunctioning equipment should no longer be used. Contact Kolibree or your retailer.
12. This product is not recommended for use by children under 3 years of age.

Medical warnings

If you are under treatment for dental problems, consult your dentist before using this product.

If you have had dental surgery in the last 2 months, consult your dentist before using this product.

Contact your dentist if you notice abnormal bleeding from gums after using this product or if bleeding persists after a week of use.

If you wear a pacemaker or an implanted device, contact your physician or the manufacturer before using this product.

Consult your doctor before using this product if you have medical problems.

Kolibree has been tested and meets the safety standards for electromagnetic devices.

12. Legal information

Modifications

The explanations and specifications contained in this user guide are provided for information purposes and may be changed without prior notification. This information is correct at the time it is published and it was written with the greatest care in order to provide you with the most accurate information possible.

Kolibree cannot be held responsible, directly or indirectly, for any damage or accidental data loss resulting from any error or omission in this document. Kolibree reserves the right to amend or improve the product and its user guide without any restriction or obligation to alert the user.

Warning and returns

Kolibree guarantees its products for one year from the date of purchase. Manufacturing defects will be repaired or replaced by Kolibree if proof of purchase can be provided and if defects are found during the warranty period.

Disclaimer of warranty:

The warranty does not cover brushes, damage caused by unintended use, user negligence, or unauthorized repairs.

Contact our customer service for more information.

Important product information:

For safety, compliance, recycling, and other important information regarding your Kolibree toothbrush, please refer to the Kolibree toothbrush user's guide which is accessible on your Kolibree application.

Kolibree: License Agreement And Terms Of Use

BEFORE USING THE KOLIBREE TOOTHBRUSH DEVICE, PLEASE READ THE LICENSE AGREEMENT AND TERMS OF USE FOUND AT WWW.KOLIBREE.COM AND ALL RULES AND POLICIES FOR THE KOLIBREE TOOTHBRUSH DEVICE AND SERVICES RELATED TO THE DEVICE (INCLUDING, BUT NOT LIMITED TO, THE KOLIBREE PRIVACY NOTICE AND ANY RULES OR USAGE PROVISIONS ACCESSIBLE THROUGH THE WWW.KOLIBREE.COM WEBSITE OR THE KOLIBREE STORE) (COLLECTIVELY, THIS "AGREEMENT"). BY USING THE KOLIBREE TOOTHBRUSH DEVICE YOU AGREE TO BE BOUND BY THE TERMS OF THIS AGREEMENT.

Contact Information. For communications concerning this Agreement, you may contact Kolibree by writing to Kolibree SAS, 36 rue de Turin, 75008 Paris, France.

Limited warranty:

This Warranty for the Kolibree toothbrush Device (the "Device") is provided by Kolibree SAS, 36 rue de Turin, 75008 Paris, France.

We warrant the Device against defects in materials and workmanship under ordinary consumer use for a period of one year or the current legal period after the initial date and place of purchase. During this warranty period, if a defect arises in the Device, and you follow the instructions for returning the Device, we will, at our option and to the extent permitted by law, either (i) repair the Device using either new or refurbished parts, (ii) replace the Device with a new or refurbished Device, or (iii) refund to you all or part of the purchase price of the Device. This limited warranty applies to any repair, replacement part or replacement Device for the remainder of the original warranty period or for ninety days, whichever period is longer. All replaced parts and Devices for which a refund is given shall become our property. This limited warranty applies only to hardware components of the Device that are not subject to accident, misuse, neglect, fire or other external causes, alterations, repair, or commercial use.

Instructions. For specific instructions about how to obtain warranty service for your Device, please contact Customer Service using the contact information provided in this guide. In general, you will need to deliver your Device in either its original packaging or in equally protective packaging to the address specified by Customer Service. Before you deliver your Device for warranty service, it is your responsibility to back up any data, software, or other materials you may have stored or preserved on your Device. It is possible that such data, software or other materials will be lost or reformatted during service, and we will not be responsible for any such damage or loss.

Limitations. TO THE EXTENT PERMITTED BY LAW, THE WARRANTY AND REMEDIES SET FORTH ABOVE ARE EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES AND REMEDIES, AND WE SPECIFICALLY DISCLAIM ALL STATUTORY OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND AGAINST HIDDEN OR LATENT DEFECTS. IF

WE CANNOT LAWFULLY DISCLAIM STATUTORY OR IMPLIED WARRANTIES, THEN TO THE EXTENT PERMITTED BY LAW, ALL SUCH WARRANTIES SHALL BE LIMITED IN DURATION TO THE DURATION OF THIS EXPRESS LIMITED WARRANTY AND TO REPAIR OR REPLACEMENT SERVICE.

SOME JURISDICTIONS DO NOT ALLOW LIMITATIONS ON HOW LONG A STATUTORY OR IMPLIED WARRANTY LASTS, SO THE ABOVE LIMITATION MAY NOT APPLY TO YOU. WE ARE NOT RESPONSIBLE FOR DIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY BREACH OF WARRANTY OR UNDER ANY OTHER LEGAL THEORY. IN SOME JURISDICTIONS THE FOREGOING LIMITATION DOES NOT APPLY TO DEATH OR PERSONAL INJURY CLAIMS, OR ANY STATUTORY LIABILITY FOR INTENTIONAL AND GROSS NEGLIGENT ACTS AND/OR OMISSIONS, SO THE ABOVE EXCLUSION OR LIMITATION MAY NOT APPLY TO YOU. SOME JURISDICTIONS DO NOT ALLOW THE EXCLUSION OR LIMITATION OF DIRECT, INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE EXCLUSION OR LIMITATION MAY NOT APPLY TO YOU. THIS "LIMITATIONS" SECTION DOES NOT APPLY TO CUSTOMERS IN THE UNITED KINGDOM AND THE EUROPEAN UNION.

This limited warranty gives you specific legal rights. You may have additional rights which vary from jurisdiction to jurisdiction, and this limited warranty does not affect those rights.


You can see a version of the Kolibree License Agreement and Terms of Use, One-Year Limited Warranty and other Kolibree information in other languages at <http://www.kolibree.com>.

13. Sending comments

Sending comments

If you want to send a message to the Kolibree team to give your opinion or ask a question, we invite you to do so by clicking on our website's ["Contact"](#) tab.

From the mobile app, if you want to send us feedback or report a bug, shake your device and follow the instructions that appear.

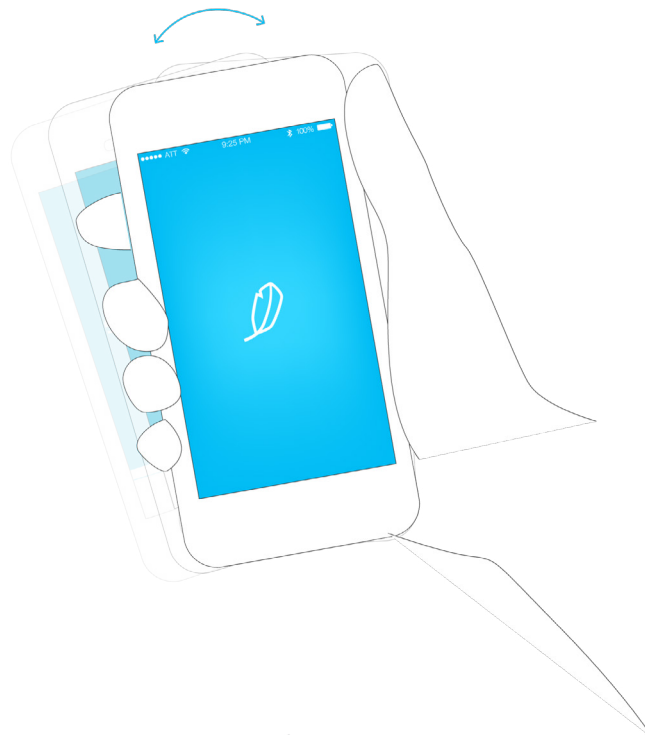

Designed by Kolibree in France. Assembled in the Philippines.

Kolibree © 2016 All rights reserved.

iPhone, iPod and iPad are registered trademarks of Apple Inc.
registered in the US and other countries. Android is a trademark of
Google Inc.